

Direct and Simultaneous Measurement of the Stiffness and Internal Friction of a Single Folded Protein


Surya Pratap S. Deopa, Shatruhan Singh Rajput, Aadarsh Kumar, Shivprasad Patil*

Department of Physics
Indian Institute of Science Education & Research
Pune 411008, Maharashtra, India
s.patil@iiserpune.ac.in

ABSTRACT

The nanomechanical response of a folded single protein, the natural nanomachine responsible for myriad biological processes, provides insight into its function. The conformational flexibility of a folded state, characterized by its viscoelasticity, allows proteins to adopt different shapes to perform their function. Despite efforts, its direct measurement has not been possible so far. We present a direct and simultaneous measurement of the stiffness and internal friction of the folded domains of the protein titin using a special interferometer based atomic force microscope. We analysed the data by carefully separating different contributions affecting the response of the experimental probe to obtain the folded state's viscoelasticity. Above ~ 95 pN of force, the individual immunoglobulins of titin transition from an elastic solid-like native state to a soft viscoelastic intermediate.

TOC figure

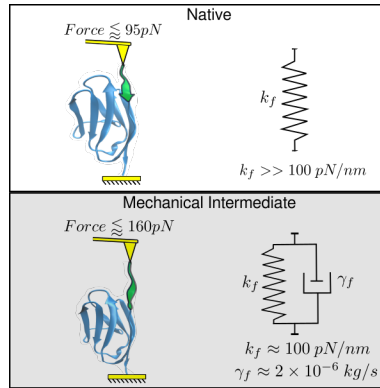

Proteins are a type of heteropolymers, which play a pivotal role in biological processes responsible for life.¹⁻³ They are natural nanomachines¹ evolved to achieve specific tasks in cellular processes.^{2,4} Proteins need to be dynamic⁵ and flexible.⁶ They also exhibit heterogeneity revealed by single molecule experiments.^{7,8} One of the important goals of nanotechnology research is to quantify folded protein's mechanical response, establish its connection with constituent network of bonds and use this knowledge for design and fabrication of efficient nanomachines.^{9,10} Owing to heterogeneity in folded states, it is important that nanomechanical response is measured at the level of single proteins.⁷

From a nano-mechanics standpoint, internal friction and stiffness can be directly measured by deforming a single protein at a certain strain rate. To characterize the viscoelasticity of a single protein one needs to strain it periodically and measure the in-phase and out-of-phase components of the stress generated. Such measurements provide both, stiffness and internal friction coefficient of the molecule. Although viscoelastic behaviour in an ensemble of proteins

has been reported,¹¹ we lack a quantitative nanometrology tool that can simultaneously measure stiffness and internal friction and is able to provide full characterization of viscoelasticity of a single protein's folded state.

It is possible to deform a single folded protein and measure the resulting stress using AFM.⁶ The conventional protein pulling experiments measure force required to unfold a single domain out of multiple domains covalently connected to each other.¹² In a typical saw-tooth pattern of such a force-extension curve, the profile between two peaks is that of the entropic elasticity of unfolded chain which is around ~ 10 pN/nm.^{13,14} The folded domains are relatively stiffer (~ 200 pN/nm)^{5,6} and hence are considered not to contribute to measured response. At a high stretch, where the folded domains stiffness is comparable to the unfolded chain, one is still faced with the challenge of separating them to obtain the response of the folded domains alone. It has been pointed out that pulling experiments using AFM are not suitable to measure viscoelasticity of the folded state.¹⁵

In this letter, we report quantitative measurement of viscoelasticity of the folded state of a single protein using a novel interferometer based AFM. By careful analysis, we have separated the response of individual folded domains from other contributions. The experiments demonstrate for the very first time that direct and simultaneous measurement of internal friction of folded protein is possible. The salient features of our instrument are i) it measures displacement of the cantilever and not bending ii) allows use of stiff cantilevers (≈ 0.6 N/m), small amplitudes (0.1 to 1 nm) and iii) in off-resonance condition, performs artefact-free measurement of nano-scale viscoelasticity in liquid environments. See Fig. 1. Briefly, a fiber with semi-mirror at its end is aligned perpendicular to a cantilever using a five-axis fibre nanopositioner (Fig.1 (a)). The mirror and back of the cantilever form a Fabry-Perot etalon. Infrared light(1310 nm) from a laser diode is guided into a 2×2 splitter. 50 % laser light goes to the etalon. The resulting interference signal, due to primary beam entering into the etalon and its interference with multiple reflections between mirror and cantilever's back surface is guided using the same splitter onto a photodiode. The current in the photodiode is extremely sensitive to fibre cantilever distance. We typically obtain sensitivity 300-500 mV/nm and measure cantilever displacement with precision of less than an angstrom.

We used repeats of 8 domains of Immunoglobulin (IgG) of titin (I-27)₈, which unfolds sequentially. 100 μ l of I-27₈ protein solution with a concentration of 10mg/ml in PBS (pH 7.4) was drop cast onto a gold coated cover-slip mounted in a fluid cell and excess protein is washed away after incubation. The cantilever is lowered into the solution and fiber is positioned on the back of the cantilever above the tip. The semi-mirror on the fibre end is aligned parallel to the back-surface of cantilever. The cantilever base is oscillated with amplitude A_0 at a frequency ω , which is far below its resonance. The amplitude A and phase lag θ of the protein-attached tip is recorded using the interferometer, as the mean separation between the base and substrate(z) is varied. The amplitude, A and phase θ versus mean separation z is used to calculate the stiffness and friction coefficient profiles shown in Fig. 2, using equations $\bar{k} = k_c \left(\frac{A_0}{|A|} \cos \theta - 1 \right)$ and $\bar{\gamma} = \frac{k_c A_0}{|A| \omega} \sin \theta$. See supplementary materials, section 1 and 2. At relatively large amplitudes ($A_0 \sim 1$ nm), and high stretch (above ~ 95 pN), a portion of the tip amplitude (~ 0.1 nm) is contributed by the folded domains.

Figure 2 shows measurement of stiffness (continuous blue line) and friction coefficient (continuous red line) of the polypeptide (I-27)₈ as the separation z between the substrate and the cantilever base is increased in a quasi-static manner. It is observed that for the first four unfolding domains, after the initial increase at lower extensions, the stiffness abruptly goes down or plateaus and rises again(the shaded region in fig. 2). This is not seen in the last four domains. The total measured dissipation $\bar{\gamma}$ also decreases with each consecutive unfolding peak. Both observations are explained in Supplementary Information, section 3.)

The Marko-Siggia approximation describes the force extension behaviour of unfolded domains in conventional experiments.¹⁶ It provides a force-extension relationship, which can be fitted to experimental force curves to obtain relevant parameters such as persistence length p and contour length L_c . The derivative of force with respect to extension provides a relationship between local stiffness k_{wlc} and end-to-end molecular extension z .¹⁶

$$k_{wlc} = \frac{dF}{dz} = \frac{k_b T}{p L_c} \left(\frac{1}{2(1 - z/L_c)^3} + 1 \right) \quad (1)$$

At lower extensions in each unfolding profile of first four domains, where stiffness of chain is much smaller compared to folded domains, wlc model of entropic elasticity given by equation 1 fits well (green continuous line). In the shaded region, however, the measured stiffness (continuous blue line) is lower than entropic stiffness. Corresponding to this region, phase also shows lag and there is measurable dissipation of energy (continuous red line).

We argue that in the shaded region of Fig.2, measured stiffness contains contribution from folded domains. We developed a method to separate it from measured viscoelastic response in this region (See supplementary information, section 1 and 2.) The conventional constant velocity pulling experiments using AFM have measured mechanical stability of I-27 domains. It mechanically unfolds via an intermediate. The domains fall into mechanical intermediate above ~ 95 pN and completely unfold at ~ 200 pN.¹⁷ See Fig. 3. Note that when the response deviates from wlc description

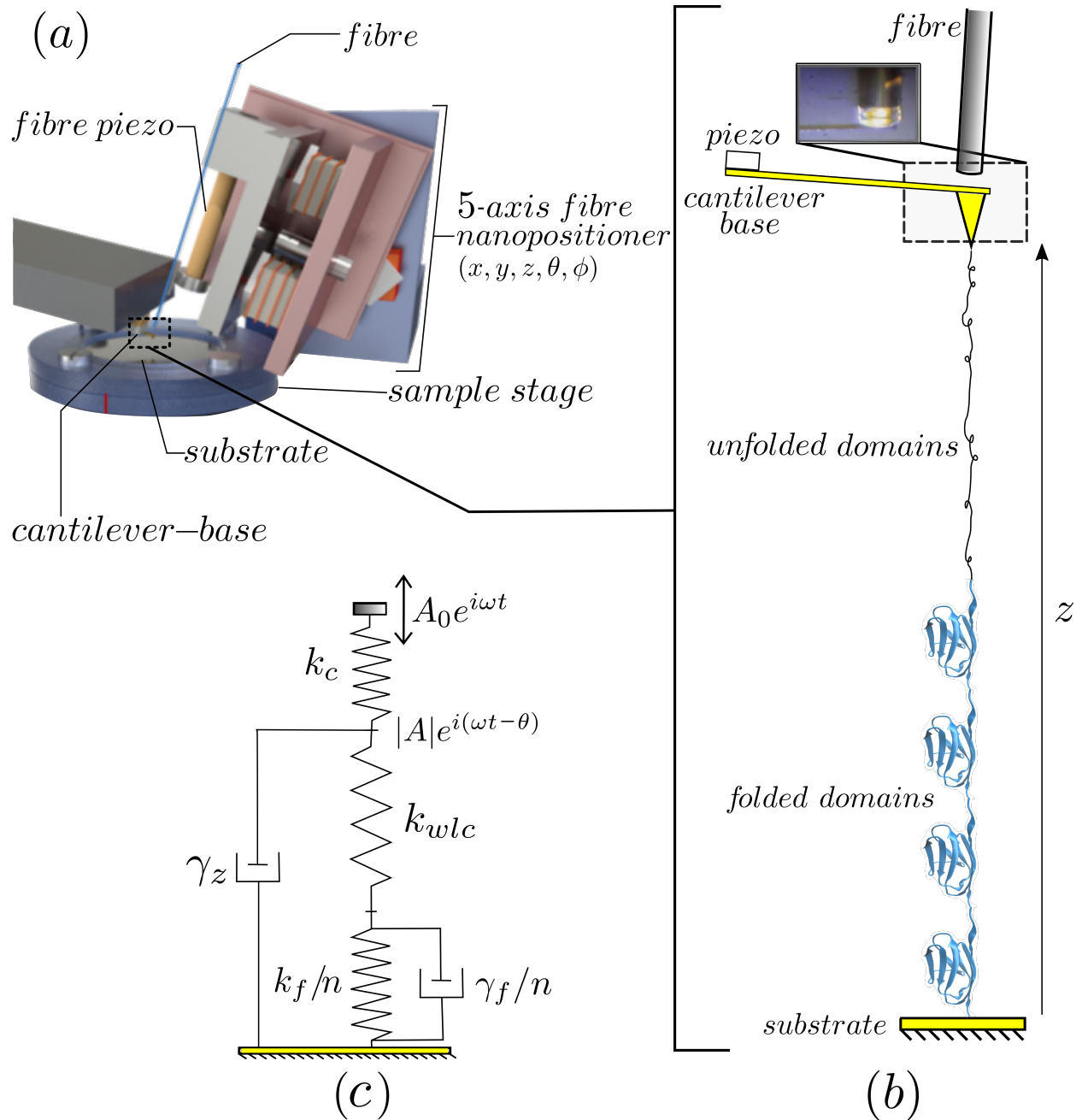


Figure 1: Measurement of folded domain's viscoelasticity using interferometer-based AFM. a) Schematic of the nano-positioner used to place the fibre over the end of the cantilever beam. The end of the fibre is coated with a semi-mirror and it is aligned parallel to the back surface of the cantilever using a five-axis inertial slider. The amplitude of the molecular extension is directly measured as opposed to the conventional method in which cantilever bending is measured. b) The measurement scheme depicting the octamer tethered to the cantilever and substrate. The base is oscillated at off-resonance while the substrate is pulled away in a quasi-static manner to increase tip-sample separation z and the amplitude (A) and phase lag (θ) are monitored. The inset shows the cleaved end of fiber coated with a semi-mirror aligned perpendicular to the length of the cantilever c) The cantilever-protein system is modelled with spring and dash-pots. A_0 - base amplitude, A - tip amplitude, ω - drive frequency, θ - phase lag between tip and base of the cantilever, k_c - cantilever stiffness γ_z hydrodynamic damping due to the proximity of the cantilever to the substrate, k_f and γ_f -stiffness and friction coefficient of folded domain respectively, n - number of folded domains. k_{wlc} - stiffness of the unfolded Worm-Like-Chain (wlc). Since k_f/n and k_{wlc} are in series, the weaker among the two springs has a dominant contribution to the amplitude A .

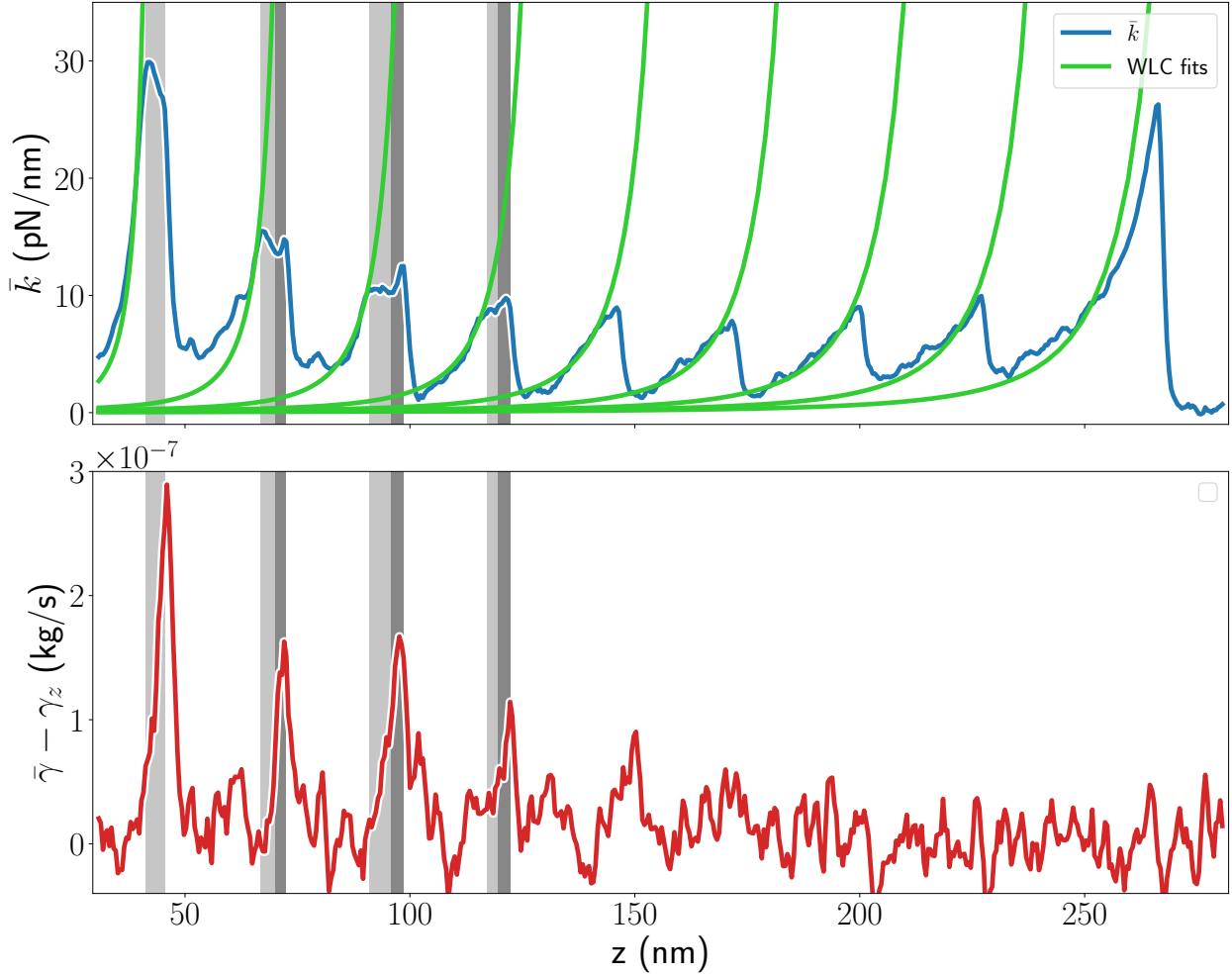


Figure 2: The measured combined stiffness (blue continuous) and friction coefficient (red continuous) of unfolded chain and the folded domains of I-27 octamer as domains sequentially unfold. The amplitude and phase of the tip displacement is used for calculation. The data shows a sawtooth pattern of unfolding events similar to constant velocity pulling experiments. The green continuous line is a fit using wlc model. The mean difference in contour length between two consecutive peaks is 29 ± 0.8 nm. A persistence length of 0.4 nm is used. The first four events deviate from wlc in the shaded region where folded domains are comparable in stiffness to the unfolded chain and contribute to the measurement. The dissipation is seen in the corresponding region due to the internal friction in folded domains. The data in the dark shaded region is used for further analysis to obtain the stiffness and internal friction coefficient of the folded state. Cantilever stiffness is 0.6 N/m, The cantilever base is excited with frequency 2.1 kHz and amplitude 1 nm. The plot is representative of more than 50 traces recorded at different frequencies and amplitudes. The procedure to obtain the molecular friction coefficient by separating the contribution from the cantilever damping at close proximity with the substrate is provided in supplementary materials, section 2.

in Fig.2, a folded domain is likely to be in either native or intermediate state and hence number of domains in either state is not known. However, as z is increased, loading the folded domains further, all of them fall into intermediate state, each adding 0.66 nm of contour length to the chain.¹⁷⁻¹⁹ At this point, as seen in figure 2, the stiffness starts to rise again before the next domain is unfolded. In this region, shaded in dark, we use our analysis method to obtain the stiffness and friction coefficient of the folded domains.

In our experiments, the oscillating base of the cantilever and the substrate are pulled away from each other with protein octamer I-27₈ tethered between the cantilever tip and the substrate (Fig.1). The measured response of the tip (amplitude A and phase θ) is due to i)that of the individual folded domains, ii) worm-like chain of unfolded domains, iii) freely jointed chain formed of folded domains. However, the freely jointed chain is 95% extended and is much stiffer compared to the first two at forces above 20 pN. Hence it does not contribute to the modulation of tip amplitude

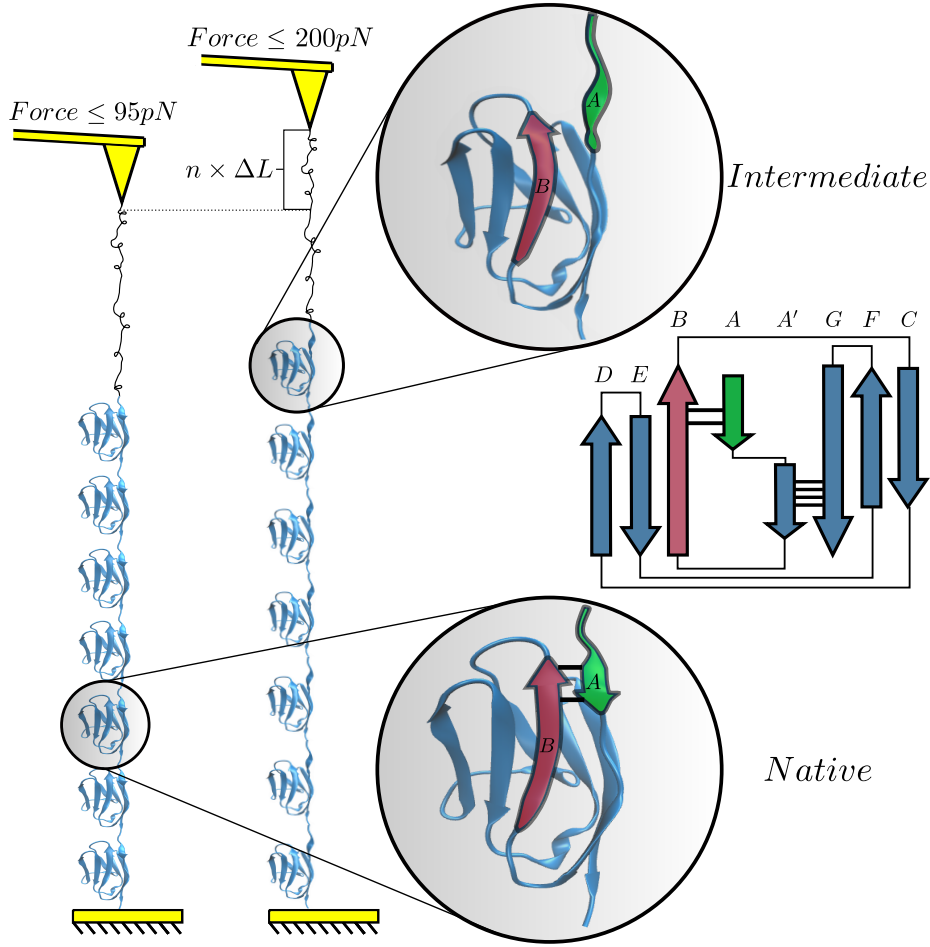


Figure 3: The octomer of I-27 under force. The cartoon shows a situation, wherein one domain is unfolded and the remaining seven are in the folded state. When the folded domains are pulled on by increasing z , the force is applied to domains through wlc chain of first unfolded domain. When force exceeds ≈ 95 pN (shaded region in fig. 2), the domains in the native state (left) are pushed into mechanical intermediate(right). As shown in the β -sheet schematic, the hydrogen bond network between B-A and A'-G is responsible for the mechanical stability of the native state, of which the network between A and B is broken for the intermediate, elongating the protein by $n \times \Delta L$, where $\Delta L = 0.66$ nm.¹⁷ When all the domains enter into the intermediate, the stiffness starts to rise again as shown in the darker shaded area in figure 2. Measurements in this region are used to estimate the stiffness and friction coefficient of the individual folded intermediates.

and only adds a constant length to the unfolded chain. We need to separate the contribution from wlc chain of unfolded domains to the measured stiffness in order to single out the viscoelastic response of folded domains. We discuss the analysis procedure followed to obtain stiffness and friction coefficient of individual folded domains. For a given z , the force balancing equation for the spring and dash-pot arrangement of Fig.1 gives us

$$k_c(A_0 - A) = \left[\left(\frac{n}{k_f + i\gamma_f\omega} + \frac{1}{k_{wlc}} \right)^{-1} + i\gamma_z\omega \right] A \quad (2)$$

Where k_{wlc} is chain stiffness and k_f and γ_f are stiffness and friction coefficient of the folded domain respectively. γ_z is damping provided by the surface to the cantilever body. We get expression for stiffness and friction coefficient of individual folded monomers as (See Supplementary Information, section 1 and 2 for details of this calculation)

$$k_f = \frac{n \left(\frac{\bar{k}}{k^2 + \omega^2\gamma^2} - \frac{1}{k_{wlc}} \right)}{\left(\frac{\bar{k}}{k^2 + \omega^2\gamma^2} - \frac{1}{k_{wlc}} \right)^2 + \left(\frac{\gamma\omega}{k^2 + \omega^2\gamma^2} \right)^2} \quad (3)$$

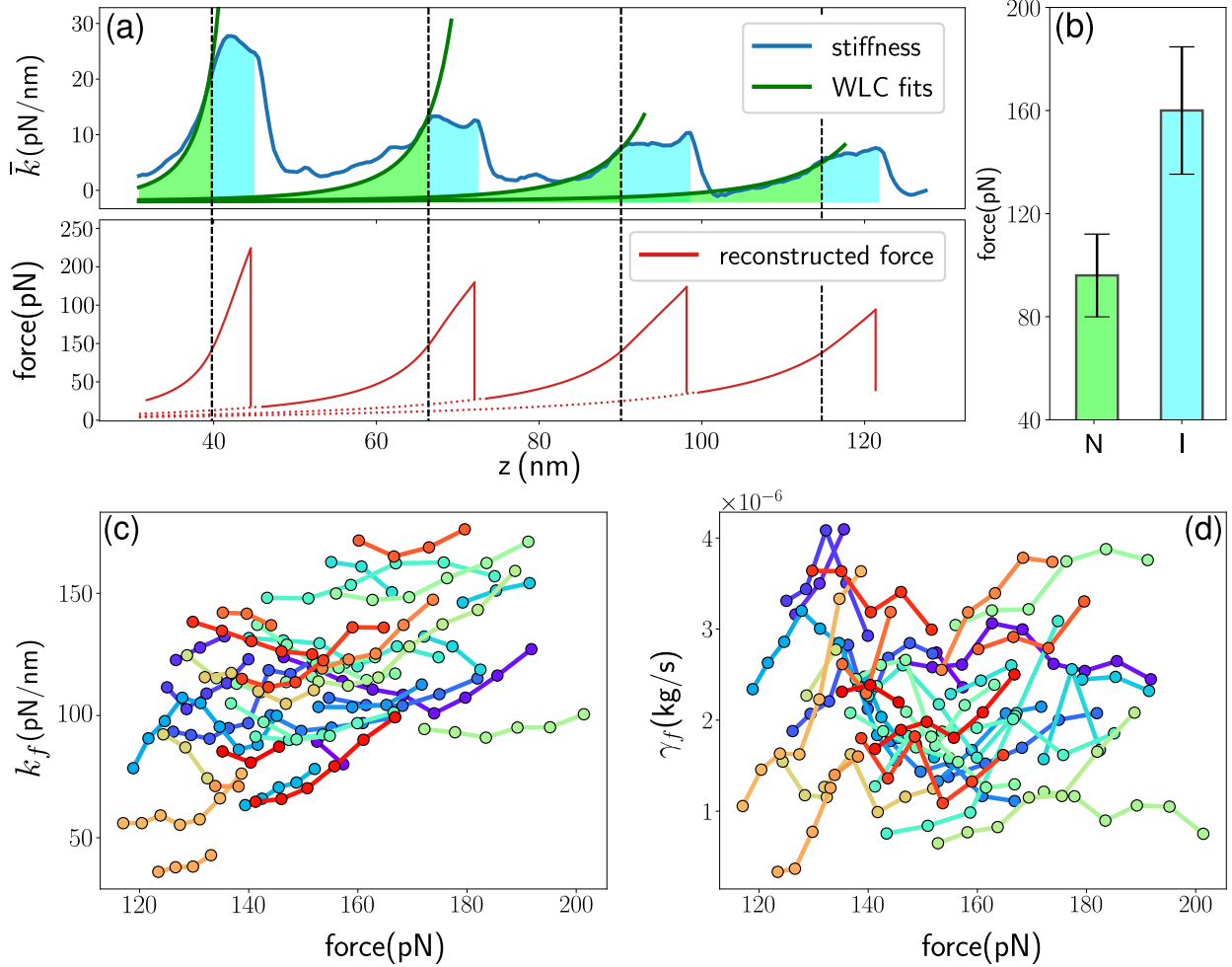


Figure 4: a) The integration of stiffness-extension (continuous blue line) provides force-extension plot (continuous red line). The vertical dotted line indicates deviation of stiffness-extension from wlc model fit (continuous green lines) of equation 1. Its intersection with force-extension gives force at which, stiffness has a contribution from folded domains. b) Bar plots of force where transition occurs and the force required to unfold a domain completely (N: Native and I: Intermediate). The deviation kicks-in when force exceeds 95 ± 16 pN. The force required to unfold the domain completely is 161 ± 25 pN. These values match with force needed to push domains of $(I-27)_8$ into the intermediate and then unfold it completely.¹⁷ The mean force required to unfold a domain completely, is ~ 200 pN in static pulling experiments. It is ~ 160 pN in our experiments due to the energy provided to it through oscillations. c) stiffness k_f and friction coefficient γ_f of a single folded domain at different forces F . The plots provide viscoelastic characterization of single folded intermediate. The data points in the same colour are from the analysis of same unfolding peak. The data is pulled together from many different profiles similar to the one shown in Fig. 2.

$$\gamma_f = \frac{n \left(\frac{\gamma \omega}{k^2 + \omega^2 \gamma^2} \right)}{\left(\frac{\bar{k}}{k^2 + \omega^2 \gamma^2} - \frac{1}{k_{wlc}} \right)^2 + \left(\frac{\gamma \omega}{k^2 + \omega^2 \gamma^2} \right)^2} \quad (4)$$

where \bar{k} and $\bar{\gamma}$ are measured stiffness and friction coefficient, $\gamma = \bar{\gamma} - \gamma_z$. We need an estimate of k_{wlc} to use equations 3 and 4 in order to get viscoelasticity of a folded domain from measured \bar{k} and $\bar{\gamma}$. For this, we integrated the measured " $\bar{k} - z$ " curve in Fig. 2 and obtained "Force- z " curve (Fig. 3). In the grey shaded region of Fig. 2, this force produces extensions both in wlc chain and the folded domains. The $k_{wlc} - F$ relationship in high force limit ($F \gg \frac{k_B T}{p}$) where the polymer is significantly extended is given by¹⁴ - $k_{wlc}(F) = \frac{4}{L_c} \sqrt{\frac{p}{k_B T}} F^{\frac{3}{2}}$. The k_{wlc} is then estimated since the

force on the chain for each z is determined. Equations 3 and 4 are used to calculate stiffness and friction coefficient of individual intermediates from experimental data.

The numerical integration of stiffness-extension data, presented in Fig.2, provides force-extension curve as shown Fig.4a. We use only first four domains from Fig. 2 for this analysis. Such force-extension profiles inform forces at which folded domains start to contribute to stiffness measurement. This is shown by dotted vertical lines intersecting the red curve depicting the force profile. The force at this intersection is threshold force after which, folded domain's stiffness contributes to measurement. This occurs at 95 ± 16 pN. The protein unfolds completely at 161 ± 25 pN. Fig. 4 b shows bar plots obtained by this procedure for all our experimental data. Using dark shaded region in Fig. 2 for our analysis, we computed stiffness (k_f) and friction coefficient (γ_f) of individual folded domains at different forces shown in Fig. 4c and d.

Fig. 4b indicates that folded protein's response starts to contribute to stiffness measurement when force on it exceeds ~ 95 pN. This matches well with force required to make the transition from native to mechanical intermediate.¹⁷⁻¹⁹ It suggests that our measurement of viscoelasticity is that of the folded intermediate, which is soft. The domains remain in the native state below ≈ 95 pN. In this region, which is not shaded in gray in Fig. 2a, the data fits wlc behavior. It indicates that in this region folded domains in the native state are too stiff to make any discernible effect on tip amplitude. It is fully elastic having immeasurably high stiffness compared to the unfolded chain. Under the loading force F , the mechanical intermediate's stiffness increases slightly with force, whereas the friction coefficient does not show any appreciable change.

It has not been possible to measure viscoelastic response of single folded protein using AFM, a technique at the forefront of pulling single proteins. The conventional deflection detection scheme, used to measure cantilever response, records cantilever bending as opposed to displacement in our experiments. (See Supplementary Information, section 4). It is important to operate off-resonance in order to avoid the stiffness change in the molecule affecting the phase, which produces measurement artefacts. The bending signal measured with deflection detection scheme is immeasurably low in off-resonance operation (See Supplementary Information, Fig. S12.)

It is possible to estimate Young's modulus G and the viscosity η from the measurement of stiffness and friction coefficient using dimensional analysis argument. $G = k \times L/A$, Where L is length, A is the cross sectional area of the object and k is its stiffness. Since we do not know the orientation of protein, taking it to be a cube of side 4 nm, results in an estimate of $G \sim 0.2$ GPa. Note that the native state has an even higher value of G . G_{glass} is ~ 50 GPa and G_{steel} is ~ 200 GPa. For the estimate of η , we take an object of \sqrt{bd} moving in a medium of viscosity η . d is extension in the folded domain. $\gamma = 6\pi\eta\sqrt{bd}$; $\eta \sim 150$ Pa.s. This is in the same order with earlier estimates by group of Hansma from tip relaxation on Lysozyme,²⁰ however a factor of $\sim 10^5$ times higher than the ones estimated from folding-refolding experiments,²¹ which Wang and Zocchi attributed to shear thinning, $\eta \sim 1/\omega^2$.¹⁵ It is likely that, using the folding rate dependence on solvent viscosity, one measures internal friction of the transition state, rather than that of the folded structures, as in this work.

Considering the folded intermediate as a viscoelastic element, its relaxation time is $\tau = \gamma_f/k_f \approx 10 \mu s$. The dissipation of energy by intermediate at μs timescales suggests that the oscillatory perturbation provided to reaction co-ordinate couples with a slowly relaxing mode. The all-atom simulations have shown that structural relaxation processes involving the protein backbone are in the microsecond time scale.²² The high stiffness and lack of dissipation of the native state suggest that, in this case, modes coupled to the reaction coordinate relax much faster. This drastic change in dynamics by merely breaking the hydrogen bond network between strands A and B, keeping the structure otherwise intact, is an intriguing new aspect of protein dynamics.

It has been shown in pulling experiments that transition from native to mechanical intermediate of I-27 does not depend on temperature or pulling speed, whereas the unfolding force of the intermediate shows this dependence clearly.^{18,19} Based on this, it has been argued that energy landscape of the intermediate is rugged whereas the native state of I-27 works as a force buffer to protect it from large physiological forces¹⁸. Our measurements are consistent with these findings. We do not see any dissipation in the folded native state and stiffness drops from an immeasurably high value to ~ 100 pN/nm as the protein transitions from native to mechanical intermediate accompanied by dissipation, indicating ruggedness.

Borrowed from the theory describing amorphous glassy systems,²³ the energy landscape description used in the study of protein's folding-refolding dynamics is empirically successful, at least for small proteins. However, this community debates the question - to what extent the projection of multi-dimensional landscape onto a single reaction co-ordinate is accurate.²⁴ The success of one dimensional picture relies on the separation of fluctuation time-scales of the chosen co-ordinate with the rest.²⁴ We have provided the first direct observation of internal friction in a folded state through stress-strain analysis. Our measurement scheme provides a direct test of such coupling between the end-to-end distance co-ordinate and other dimensions in the reduced, multidimensional protein landscape. The internal friction measured

here is a result of energy loss in driving the folded state through many of its substates and it may provide useful insights into the inter-conversion rates.²⁵ The observed dissipation in our experiments is not possible unless the end-to-end distance reaction coordinate, which is modulated by the tip at ~ 2 KHz is coupled to slowly relaxing backbone dynamics.²² The forced oscillations of end-to-end distance at $\sim ms$ timescale results in dissipation through such coupling. The data on lower operational frequency is presented in Supplementary Information, section 5.

The protein of choice in our experiments is not only a paradigm in pulling experiments, but also provides us a native elastic state and a viscoelastic mechanical intermediate in one experiment. We show here that only the network of hydrogen bonds between β strands A and B decide the fate of mechanics governing IgG domains. The domain behaves like a Hookean solid, having no dissipation with the network intact, whereas it has viscous elements, if the network is broken. Understanding effects of such local changes on global properties of folded macromolecules are crucial for unraveling the mechanism of allostery^{26,27} and could provide an instructive design principle for synthesizing nanomachines using molecular dynamics simulations.⁹

While our stiffness measurements match well with previous efforts,^{5,6} the friction coefficient of single protein's folded state are not yet measured. We discuss new possibilities opened with this ability. Internal friction implies a rough energy landscape and slow dynamics due to frustration. Our direct measurement of dissipation also reveals that intermediate state is frustrated. It allows estimate of effective diffusion coefficient through Stoke-Einstein relationship $D = K_B T / \gamma$, where γ is the friction coefficient. From our measurement of γ for the intermediate, the diffusion coefficient is $D = (1.94 \pm 0.05) \times 10^3 \text{ nm}^2/\text{s}$. It is possible to estimate ruggedness (ϵ) by using Zwanzig relationship and measuring D at two different temperatures.²⁸ Secondly, using our experimental scheme, where the folded state is driven out of equilibrium, the rupture-forces distribution along with energy dissipation to maintain the non-equilibrium state can be used to extract information about the dynamics of slow processes,²⁵ which are poorly understood at present.

In conclusion, we simultaneously and directly measure stiffness and internal friction of single folded protein domain for the first time. In particular, we find that beyond a threshold force, the IgG domain of titin moves from an elastic solid-like native state to a viscoelastic mechanical intermediate, which dissipates energy, providing the first direct evidence of internal friction. The capability demonstrated in our experiment provides hope of relating mechanics to chemical details at the level of single bonds by coupling this type of measurement to suitable single-molecule optical spectroscopy, such as Förster Resonance Energy Transfer.

Supplementary Description

Model: Protein-Cantilever system; Data Analysis; Features observed in stiffness vs extension and friction vs extension profiles; Fibre Interferometer AFM and single molecule experiments; Measurement at lower frequency.

Acknowledgement

SPD acknowledges fellowship from IISER Pune. The work is supported by Wellcome-Trust DBT India alliance through intermediate fellowship to SP (500172/Z/09/Z).

References

- (1) Alberts, B.; Bray, D.; Lewis, J.; Raff, M.; Roberts, K.; Watson, J., *Molecular Biology of the Cell*, 4th; Garland, New York: 2002.
- (2) Alberts, B. The cell as a collection overview of protein machines: Preparing the next generation of molecular biologists. *cell* **1998**, *92*, 291–294.
- (3) Goodsell, D. S., *The machinery of life*; Copernicus New York, NY: 2009.
- (4) Thirumalai, D.; Hyeon, C.; Zhuravlev, P. I.; G.H., L. Symmetry, Rigidity, and Allosteric Signaling: From Monomeric Proteins to Molecular Machines. *Chemical Reviews* **2019**, *119*, 6788–6821.
- (5) Zaccai, G. How soft is a protein? A protein dynamics force constant measured by neutron scattering. *Science* **2000**, *288*, 1604–1607.
- (6) Dong, M.; Husale, S.; Sahin, O. Determination of protein structural flexibility by microsecond force spectroscopy. *Nature nanotechnology* **2009**, *4*, 514–517.
- (7) Hyeon, C.; Hinczewski, M.; Thirumalai, D. Evidence of disorder in biological molecules from single molecule pulling experiments. *Physical review letters* **2014**, *112*, 138101.
- (8) Hinczewski, M.; Hyeon, C.; Thirumalai, D. Directly measuring single-molecule heterogeneity using force spectroscopy. *Proceedings of the National Academy of Sciences* **2016**, *113*, E3852–E3861.

- (9) Soong, R. K.; Bachand, G. D.; Neves, H. P.; Olkhovets, A. G.; Craighead, H. G.; Montemagno, C. D. Powering an inorganic nanodevice with a biomolecular motor. *Science* **2000**, *290*, 1555–1558.
- (10) Howard, J., *Mechanics of motor proteins and the cytoskeleton*; Sinauer associates Sunderland, MA: 2001.
- (11) Wang, Y.; Zocchi, G. Elasticity of globular proteins measured from the ac susceptibility. *Physical review letters* **2010**, *105*, 238104.
- (12) Carrion-Vazquez, M.; Oberhauser, A. F.; Fowler, S. B.; Marszalek, P. E.; Broedel, S. E.; Clarke, J.; Fernandez, J. M. Mechanical and chemical unfolding of a single protein: a comparison. *Proceedings of the National Academy of Sciences* **1999**, *96*, 3694–3699.
- (13) Rajput, S. S.; Deopa, S. P. S.; Yadav, J.; Ahlawat, V.; Talele, S.; Patil, S. The nano-scale viscoelasticity using atomic force microscopy in liquid environment. *Nanotechnology* **2020**, *32*, 085103.
- (14) Khatri, B. S.; Byrne, K.; Kawakami, M.; Brockwell, D. J.; Smith, D. A.; Radford, S. E.; McLeish, T. C. Internal friction of single polypeptide chains at high stretch. *Faraday discussions* **2008**, *139*, 35–51.
- (15) Wang, Y.; Zocchi, G. The folded protein as a viscoelastic solid. *EPL (Europhysics Letters)* **2011**, *96*, 18003.
- (16) Marko, J. F.; Siggia, E. D. Stretching dna. *Macromolecules* **1995**, *28*, 8759–8770.
- (17) Marszalek, P. E.; Lu, H.; Li, H.; Carrion-Vazquez, M.; Oberhauser, A. F.; Schulten, K.; Fernandez, J. M. Mechanical unfolding intermediates in titin modules. *Nature* **1999**, *402*, 100–103.
- (18) Nunes, J. M.; Hensen, U.; Ge, L.; Lipinsky, M.; Helenius, J.; Grubmüller, H.; Muller, D. J. A “force buffer” protecting immunoglobulin titin. *Angewandte Chemie International Edition* **2010**, *49*, 3528–3531.
- (19) Taniguchi, Y.; Brockwell, D. J.; Kawakami, M. The effect of temperature on mechanical resistance of the native and intermediate states of I27. *Biophysical journal* **2008**, *95*, 5296–5305.
- (20) Radmacher, M.; Fritz, M.; Cleveland, J. P.; Walters, D. A.; Hansma, P. K. Imaging adhesion forces and elasticity of lysozyme adsorbed on mica with the atomic force microscope. *Langmuir* **1994**, *10*, 3809–3814.
- (21) Tseng, C.-Y.; Wang, A.; Zocchi, G.; Rolih, B.; Levine Alex J Cellmer, T.; Henry, E. R.; Hofrichter, J.; Eaton, W. A. Measuring internal friction of an ultrafast-folding protein. *Proceedings of National Academy of Sciences* **2008**, *105*, 18320.
- (22) Shaw, D. E.; Maragakis, P.; Lindorff-Larsen, K.; Piana, S.; Dror, R. O.; Eastwood, M. P.; Bank, J. A.; Jumper, J. M.; Salmon, J. K.; Shan, Y., et al. Atomic-level characterization of the structural dynamics of proteins. *Science* **2010**, *330*, 341–346.
- (23) Frauenfelder, H.; Sligar, S. G.; Wolynes, P. G. The energy landscapes and motions of proteins. *Science* **1991**, *254*, 1598–1603.
- (24) Neupane, K.; Manuel, A. P.; Woodside, M. T. Protein folding trajectories can be described quantitatively by one-dimensional diffusion over measured energy landscapes. *Nature Physics* **2016**, *12*, 700–703.
- (25) Horowitz, J. M.; Zhou, K.; England, J. L. Minimum energetic cost to maintain a target nonequilibrium state. *Physical Review E* **2017**, *95*, 042102.
- (26) Schachman, H. H. K. Still looking for ivory tower. *Annual Review of Biochemistry* **2000**, *69*, 1–69.
- (27) Weber, G. Ligand binding and internal equilibiums in proteins. *Biochemistry* **1972**, *11*, 864–878.
- (28) Zwanzig, R. Diffusion in a rough potential. *Proceedings of the National Academy of Sciences* **1988**, *85*, 2029–2030.
- (29) Rief, M.; Gautel, M.; Oesterhelt, F.; Fernandez, J. M.; Gaub, H. E. Reversible unfolding of individual titin immunoglobulin domains by AFM. *science* **1997**, *276*, 1109–1112.
- (30) Kawakami, M.; Byrne, K.; Brockwell, D. J.; Radford, S. E.; Smith, D. A. Viscoelastic study of the mechanical unfolding of a protein by AFM. *Biophysical journal* **2006**, *91*, L16–L18.
- (31) O’Shea, S. Comment on “oscillatory dissipation of a simple confined liquid”. *Physical review letters* **2006**, *97*, 179601.
- (32) Kiracofe, D.; Raman, A. Quantitative force and dissipation measurements in liquids using piezo-excited atomic force microscopy: a unifying theory. *Nanotechnology* **2011**, *22*, 485502.
- (33) Benedetti, F.; Gazizova, Y.; Kulik, A. J.; Marszalek, P. E.; Klinov, D. V.; Dietler, G.; Sekatskii, S. K. Can dissipative properties of single molecules be extracted from a force spectroscopy experiment? *Biophysical journal* **2016**, *111*, 1163–1172.
- (34) Janovjak, H.; Müller, D. J.; Humphris, A. D. Molecular force modulation spectroscopy revealing the dynamic response of single bacteriorhodopsins. *Biophysical journal* **2005**, *88*, 1423–1431.
- (35) Taniguchi, Y.; Khatri, B. S.; Brockwell, D. J.; Paci, E.; Kawakami, M. Dynamics of the coiled-coil unfolding transition of myosin rod probed by dissipation force spectrum. *Biophysical journal* **2010**, *99*, 257–262.
- (36) Smith, D. T.; Pratt, J. R.; Howard, L. A fiber-optic interferometer with subpicometer resolution for dc and low-frequency displacement measurement. *Review of Scientific Instruments* **2009**, *80*, 035105.

Supplementary Information

Contents

1 Model: Protein - Cantilever System	10
2 Data Analysis	15
2.1 Analysis to obtain k_{wlc}	15
2.2 Obtaining k_f and γ_f	16
3 Features observed in stiffness-extension and friction coefficient-extension profiles	20
3.1 The wlc plots of stiffness verses extension and Force verses extension	20
3.2 Decreasing peak in total friction coefficient for consecutive unfolding events	21
4 Fibre-Interferometer based AFM and single molecules	21
5 Measurement at lower frequency	25

1 Model: Protein - Cantilever System

There are two aspects to the experiment. Firstly, as the separation z between the cantilever base and the substrate is increased, force F builds up on the protein that is attached between the cantilever tip and the surface. Eventually, one of the protein domains unfolds due to the force. The other aspect involves dynamic measurements where the base of the cantilever is oscillated with a constant frequency and amplitude. The substrate retraction is halted at each separation z for a period of 30 ms. During this period the dynamic response of the protein under force F (due to the cantilever base and surface separation z) is obtained by recording the amplitude and phase of the cantilever tip. The equation of motion for the cantilever is given by

$$m\ddot{z} + (k_c + \bar{k})z + \bar{\gamma}\dot{z} = k_c A_0 \exp i\omega t$$

Here, z is the deflection of the cantilever, m is the effective mass of the moving cantilever, $\bar{\gamma}$ is the friction coefficient for total damping experienced by the cantilever, k_c is the cantilever stiffness, A_0 is the drive amplitude. The \bar{k} and $\bar{\gamma}$ are obtained after linearizing the conservative and dissipative components of the interaction. This assumption is valid because over the oscillation cycle of 1 nm the interaction force can be considered to be varying linearly. This is the equation of a linear forced oscillator with the well-known solutions.

$$|A| = \frac{k_c A_0}{\sqrt{(k_c + \bar{k})^2 \left(1 - \frac{\omega^2}{\omega_0^2}\right)^2 + (\bar{\gamma}\omega)^2}}$$

$$\tan \theta = -\frac{\bar{\gamma}\omega}{(k_c + \bar{k}) \left(1 - \frac{\omega^2}{\omega_0^2}\right)}$$

Where ω_0 is given by $\sqrt{\frac{k_c + k_i}{m}}$. For off-resonance regime the $1 - \frac{\omega^2}{\omega_0^2}$ term becomes unity. We can now solve for \bar{k} and $\bar{\gamma}$ to get

$$\bar{k} = k_c \left(\frac{A_0}{|A|} \cos \theta - 1 \right)$$

$$\bar{\gamma} = \frac{k_c A_0}{|A| \omega} \sin \theta$$

When operating strictly at off-resonance frequency and no molecule is attached to the tip, the phase difference between the cantilever and the drive is close to zero. The base amplitude and the tip amplitude are same. Hence we can neglect the inertial and velocity dependent forces and the response of the cantilever can be described by its static response. One can also arrive at the same expression by following a much simpler and intuitive method. For the off-resonance operation, the cantilever and the interaction are in series with each other. The conservative and dissipative components of the interaction are modelled as spring and dash-pot respectively. For such an arrangement we can write the force balance equation as

$$k_c(A_0 e^{i\omega t} - |A| e^{i(\omega t - \theta)}) = (\bar{k} + i\bar{\gamma}\omega)|A| e^{i(\omega t - \theta)} \quad (5)$$

Here A_0 and A are the base (drive) amplitude and the tip amplitude respectively. k_c is the cantilever stiffness and ω is the drive frequency. \bar{k} and $\bar{\gamma}$ are stiffness and friction coefficient of the material beneath the tip. And θ is the phase lag between the tip and the base.

$$\begin{aligned} k_c(A_0 - |A| e^{-i\theta}) &= (\bar{k} + i\bar{\gamma}\omega)|A| e^{-i\theta} \\ k_c(A_0 e^{i\theta} - |A|) &= (\bar{k} + i\bar{\gamma}\omega)|A| \\ k_c A_0 \cos \theta + i k_c \sin \theta - k_c |A| &= \bar{k} |A| + i |A| \bar{\gamma} \omega \end{aligned}$$

comparing real and imaginary components gives us

$$\bar{k} = k_c \left(\frac{A_0}{|A|} \cos \theta - 1 \right) \quad (6)$$

$$\bar{\gamma} = \frac{k_c A_0}{|A| \omega} \sin \theta \quad (7)$$

The stiffness and friction coefficient of the material can then be quantified using equations 2 and 3. To operate at off-resonance, stiff cantilevers, more than an order of magnitude higher than what are typically used in constant velocity experiments are required.

In the experiment, cantilever, folded domains and the unfolded chain all are in series as shown in Fig S1. The damping provided to the oscillating tip due to proximity of the surface is in parallel to the protein system as shown in the Fig. S1. For a particular mean tip-sample separation z , the force on all elements such as unfolded chain, folded domains and the linker connecting the domains is same.

$$k_c(A_0 - A) = \left[\left(\frac{n}{k_f + i\gamma_f \omega} + \frac{1}{k_{wlc}} \right)^{-1} + i\gamma_z \omega \right] A \quad (8)$$

Where n is number of folded domains.

Rearranging equation 4

$$\begin{aligned} k_c \left(\frac{A_0}{A} - 1 \right) &= \left(\frac{1}{k_f + i\gamma_f \omega} + \frac{1}{k_{wlc}} \right)^{-1} + i\gamma_z \omega \\ k_c \left(\frac{A_0}{|A|} e^{-i\theta} - 1 \right) &= \left(\frac{1}{k_{wlc}} + \frac{n}{k_f + i\gamma_f \omega} \right)^{-1} + i\gamma_z \omega \\ k_c \left(\frac{A_0}{|A|} \cos \theta - 1 \right) + i k_c \frac{A_0}{|A|} \sin \theta &= \left(\frac{1}{k_{wlc}} + \frac{n}{k_f + i\gamma_f \omega} \right)^{-1} + i\gamma_z \omega \end{aligned}$$

Using equations 2 and 3

$$\bar{k} + i\bar{\gamma}\omega = \left(\frac{1}{k_{wlc}} + \frac{n}{k_f + i\gamma_f \omega} \right)^{-1} + i\gamma_z \omega \quad (9)$$

$$= (\alpha - i\beta\omega)^{-1} + i\gamma_z \omega \quad (10)$$

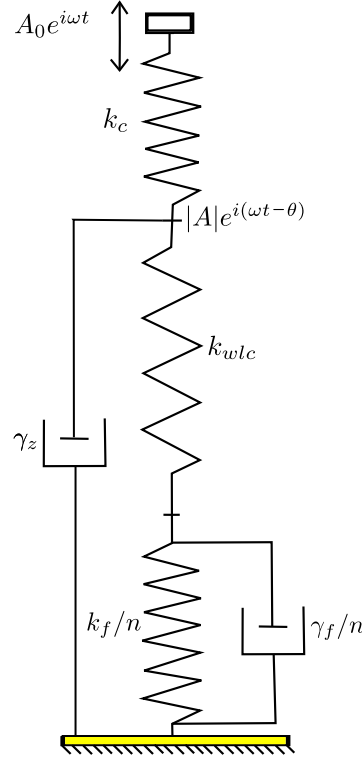


Figure S1: A spring-dashpot model for the protein-cantilever system. A_o and A are the base and tip amplitude respectively and θ is the phase lag between them. ω is the drive frequency. k_{wlc} and γ_z are the unfolded chain stiffness and surface damping respectively. k_f and γ_f are the stiffness and friction coefficient of the folded domain, where k_f/n and γ_f/n means n of them are in series.

$$\text{where } \alpha = \frac{1}{k_{wlc}} + \frac{nk_f}{k_f^2 + \omega^2\gamma_f^2} \quad \text{and} \quad \beta = \frac{n\gamma_f}{k_f^2 + \omega^2\gamma_f^2}$$

$$\bar{k} + i\bar{\gamma}\omega = \frac{\alpha}{\alpha^2 + \omega^2\beta^2} + i\omega\left(\frac{\beta}{\alpha^2 + \omega^2\beta^2} + \gamma_z\right) \quad (11)$$

The real and imaginary parts of equation can be separated to obtain measured stiffness and friction coefficient in terms of the model parameters.

$$\bar{k} = \frac{\alpha}{\alpha^2 + \omega^2\beta^2} \quad (12)$$

$$\bar{\gamma} = \frac{\beta}{\alpha^2 + \omega^2\beta^2} + \gamma_z \quad (13)$$

Case 1:

When no protein is attached to the cantilever the α and β terms in the equation 8 and equation 9 are 0. This implies that the measured stiffness, in this case, should be zero and the only damping measured should be that provided by the surface. This is seen when the cantilever is pulled off without a molecule attached to it. See figure S2(b).

$$\bar{k} = 0 \quad \bar{\gamma} = \gamma_z$$

–

Case 2:

When an octamer is picked up by the tip and all its 8 domains are pulled at, the number of folded domains $n = 8$. As domains are unfolded sequentially, "n" reduces by 1 at each unfolding event. However, if the stiffness of folded native

states is much higher than the chain stiffness ($k_{wlc} \ll k_f$), the second term in α is negligible. Similarly, β also is negligibly small.

$$\bar{k} = k_{wlc} \qquad \bar{\gamma} = \gamma_z$$

This is also supported by the experiments where the stiffness fits well to WLC model of elasticity as seen in Fig.2 in the main text.

Case 3:

It is known that in the case of I27, when it is pulled with a force above ~ 90 -100 pN, the domain makes a transition to intermediate which is softer as is evident from the data in Fig. 2 in the main text.

When n domains are in intermediate, measured stiffness and friction coefficient are given by equation 8 and equation 9 with the respective values of α and β .

$$\gamma = \bar{\gamma} - \gamma_z = \frac{\frac{n\gamma_f}{k_f^2 + \omega^2\gamma_f^2}}{\left(\frac{1}{k_{wlc}} + \frac{nk_f}{k_f^2 + \omega^2\gamma_f^2}\right)^2 + \left(\frac{n\gamma_f}{k_f^2 + \omega^2\gamma_f^2}\right)^2} \quad (14)$$

$$\bar{k} = \frac{\frac{1}{k_{wlc}} + n\frac{k_f}{k_f^2 + \omega^2\gamma_f^2}}{\left(\frac{1}{k_{wlc}} + n\frac{k_f}{k_f^2 + \omega^2\gamma_f^2}\right)^2 + \left(\frac{n\gamma_f}{k_f^2 + \omega^2\gamma_f^2}\right)^2} \quad (15)$$

Where $\gamma = \bar{\gamma} - \gamma_z$ is damping provided by the protein alone. It is obtained by subtracting out the background due to surface damping by fitting a polynomial. See Fig. S2(a).

Case 4:

With each unfolding event, the unfolded chain's contour length increases by 29 nm. The larger contour length makes the chain more compliant. The unfolding force is reached before the chain is extended to have a stiffness which is comparable to the folded intermediates (see Fig S8). In this region usually after four unfolding events, the folded intermediate's stiffness again becomes much larger compared to the unfolded chain $k_{wlc} \ll k_f$. The surface damping at such large separations can also be neglected as is seen in the data. This reduces equation 10 and 11 to

$$\bar{k} = k_{wlc} \qquad \bar{\gamma} = \gamma_z = 0$$

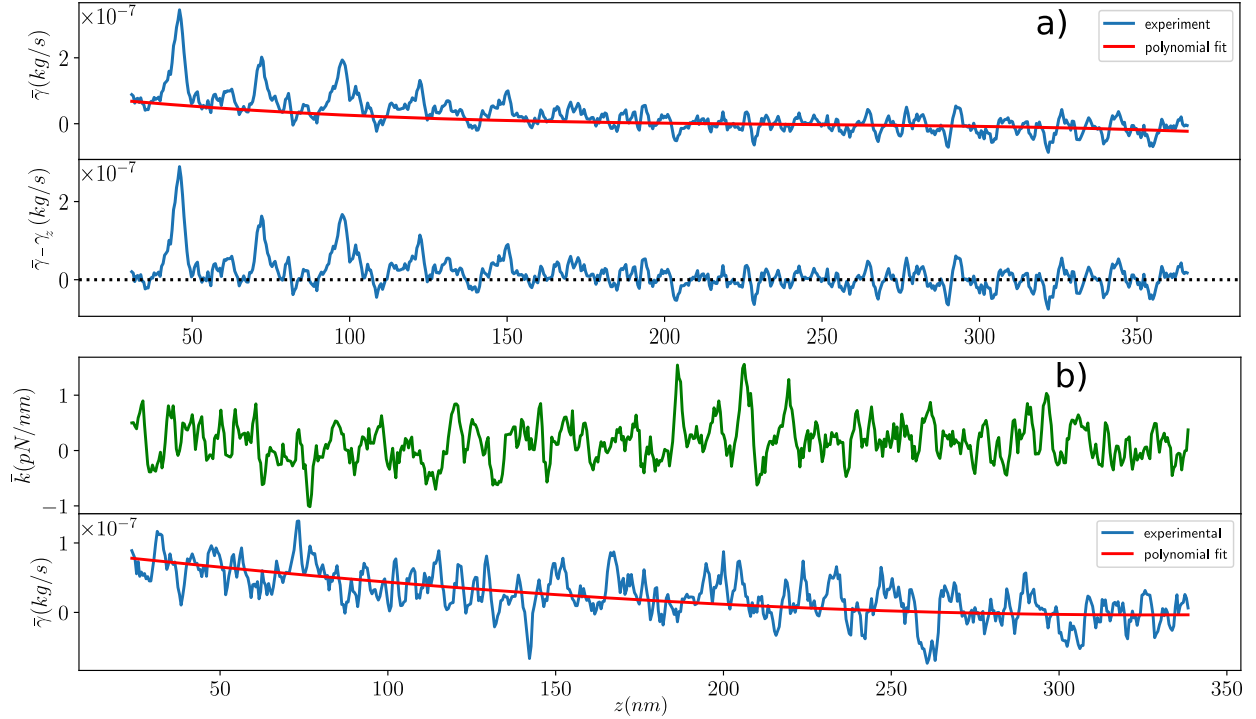


Figure S2: a) The measured total friction coefficient (blue) along with a polynomial fit (red) for the damping due to the surface effects. The lower panel in a) represents friction coefficient of the protein alone, $\bar{\gamma} - \gamma_z$ b) shows that measured stiffness is zero when the protein is not attached, however the damping due to the surface is still measured (blue).

Since our model and its special cases describe features observed in our experiments, we use equation 4 to get stiffness and friction coefficient of folded intermediates.

$$\left(\frac{1}{\bar{k} + i(\bar{\gamma} - \gamma_z)\omega} - \frac{1}{k_{wlc}} \right)^{-1} = \frac{k_f}{n} + i \frac{\gamma_f \omega}{n} \quad (16)$$

Simplifying this equation leads to the expression for stiffness and friction coefficient of the folded domains in terms of the experimentally measured quantities.

$$k_f = \frac{n \left(\frac{\bar{k}}{\bar{k}^2 + \omega^2 \gamma^2} - \frac{1}{k_{wlc}} \right)}{\left(\frac{\bar{k}}{\bar{k}^2 + \omega^2 \gamma^2} - \frac{1}{k_{wlc}} \right)^2 + \left(\frac{\gamma \omega}{\bar{k}^2 + \omega^2 \gamma^2} \right)^2} \quad (17)$$

$$\gamma_f = \frac{n \left(\frac{\gamma \omega}{\bar{k}^2 + \omega^2 \gamma^2} \right)}{\left(\frac{\bar{k}}{\bar{k}^2 + \omega^2 \gamma^2} - \frac{1}{k_{wlc}} \right)^2 + \left(\frac{\gamma \omega}{\bar{k}^2 + \omega^2 \gamma^2} \right)^2} \quad (18)$$

Where all the quantities on the right hand side are known experimentally. k_{wlc} is obtained by fitting wlc to the stiffness-extension data. In the next section, we describe our methodology to obtain k_{wlc} from our experimental data. Using equations 13 and 14, the stiffness and friction coefficient of folded intermediates are obtained when all the folded domains n are in the intermediate state. In experiments, this is marked by the dark shaded region in Fig. 2 of the main text, where the measured stiffness starts to rise again.

2 Data Analysis

2.1 Analysis to obtain k_{wlc}

Fitting the wlc model: In the initial part of rising stiffness (region in Fig. 2, not shaded, in the manuscript) for each unfolding event, all the folded domains are in the native state, where ($k_f \gg k_{wlc}$). This leads to equation 8 reducing to $\bar{k} = k_{wlc}$, as discussed in the previous section. The contribution to measured stiffness in this region is entirely from the unfolded chain, as is evident from the wlc model fits. The differential form of the wlc model¹⁶ which relates k_{wlc} to its end to end distance is fitted to the measured stiffness. The contour length L_c is the free parameter and the persistence length p is taken as 0.4 nm.^{29,30}

$$\frac{dF}{dz} = \frac{k_B T}{p L_c} \left(\frac{1}{2(1 - z/L_c)^3} + 1 \right) \quad (19)$$

It is observed that the wlc fits to the experimental data up to a certain value of the extension z , as seen in Fig. 2.

Reconstructing force profiles from stiffness data: In order to estimate the contribution from the unfolded chain, k_{wlc} , to the total measured stiffness (shaded region in Fig. 2 manuscript), we need to evaluate force in this region. This is obtained by integrating the stiffness-extension data. This integration is done in two parts. First, this integration is done over the fitted wlc curve. Second, the numerical integration is carried out on the experimental data in the region where the wlc deviates from the measured stiffness-extension profile (shaded in blue). After obtaining the force profile in the shaded blue region, one can use the following to estimate k_{wlc} . The relationship is shown to work in the region of high stretch.¹⁴

$$k_{wlc}(F) = \frac{4}{L_c} \sqrt{\frac{p}{k_B T}} F^{\frac{3}{2}} \quad (20)$$

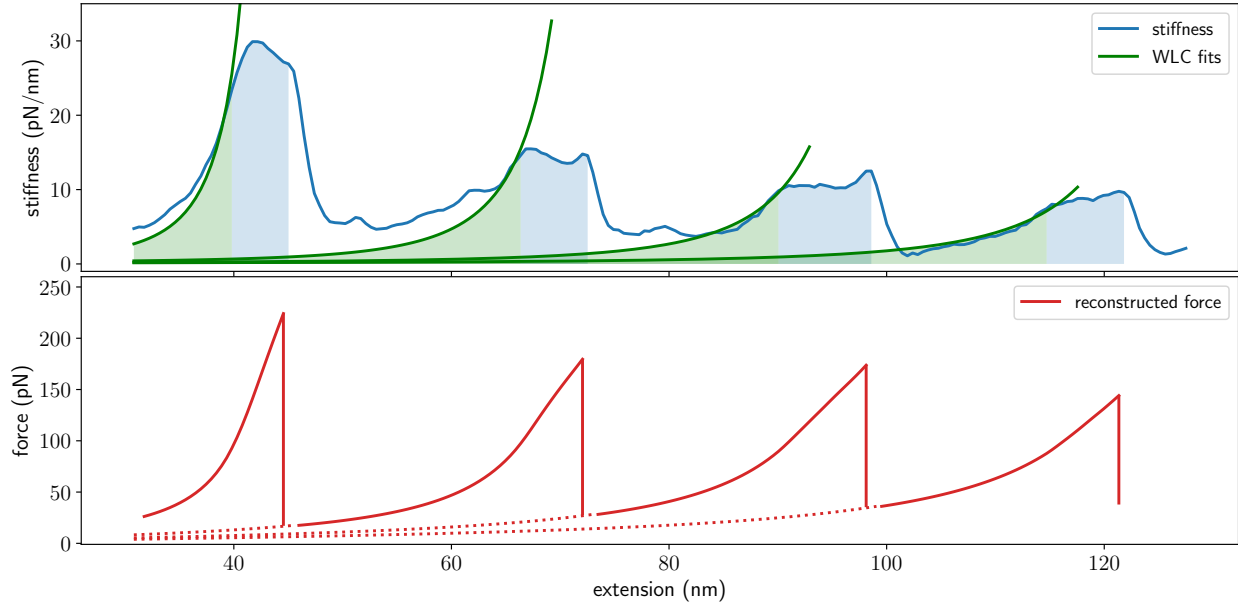


Figure S3: Numerical integration of the stiffness-extension profiles to obtain force-extension profiles. Continuous green curve- wlc fits, Continuous blue curve- experimental data, continuous red curve - force obtained after integration. The green shaded region represents integration done over fitted WLC model. The blue shaded region represents integration over experimental data where it deviates from WLC behavior.

2.2 Obtaining k_f and γ_f

In equation 16, the contour length L_c needs to be evaluated for the analysis of the shaded region. Note that in this region, some of the folded domains are in the native state and some are in intermediate. The intermediate is formed by breaking hydrogen bonds between A and B β -strands. Once broken, the β -strand A is released adding 0.66 nm to the contour length.¹⁷⁻¹⁹ As a result, it is difficult to obtain the exact value of L_c in the shaded region of a given unfolding event. Secondly, the number of domains in the intermediate n , are needed in equation 13 and 14 to obtain the stiffness of folded intermediate. One of the noteworthy features of our data is that towards the end of the shaded region the measured stiffness starts to rise again. We argue that at this point, all the folded domains are in the intermediate. This allows us to use n in equation 13 and 14 and L_c in equation 15 becomes $L_c + n0.66$. With this, the k_{wlc} at each force $F(z)$ is determined and equation 13 and 14 are then used to estimate the stiffness and friction coefficient of folded intermediate. Note that we are able to measure the stiffness of the folded intermediate using this methodology, the stiffness of the native folded state is immeasurably high.

The analysis scheme to obtain k_f and γ_f is repeated over many data sets. These data sets are shown in Fig. S4-S7. The grey shaded region is used for the said data analysis. The values obtained at different forces are then plotted together in figure 4c and d in the main text.

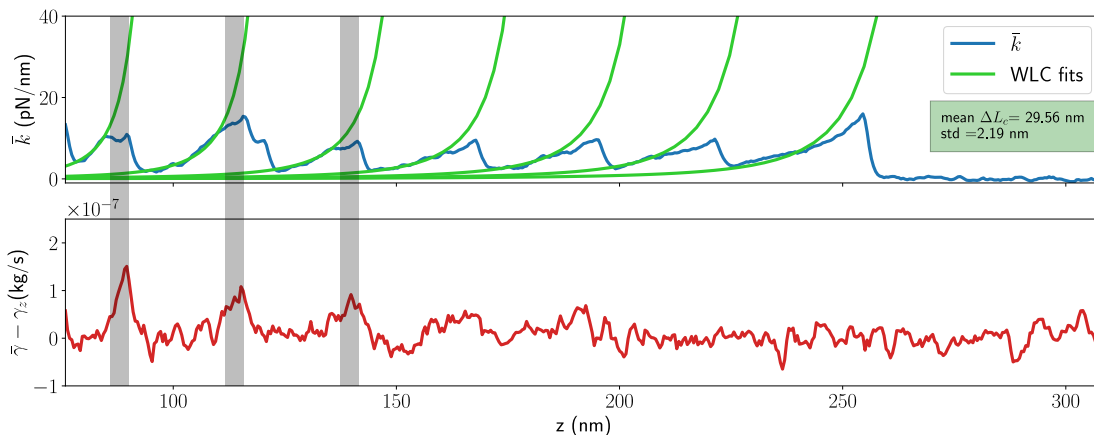


Figure S4:

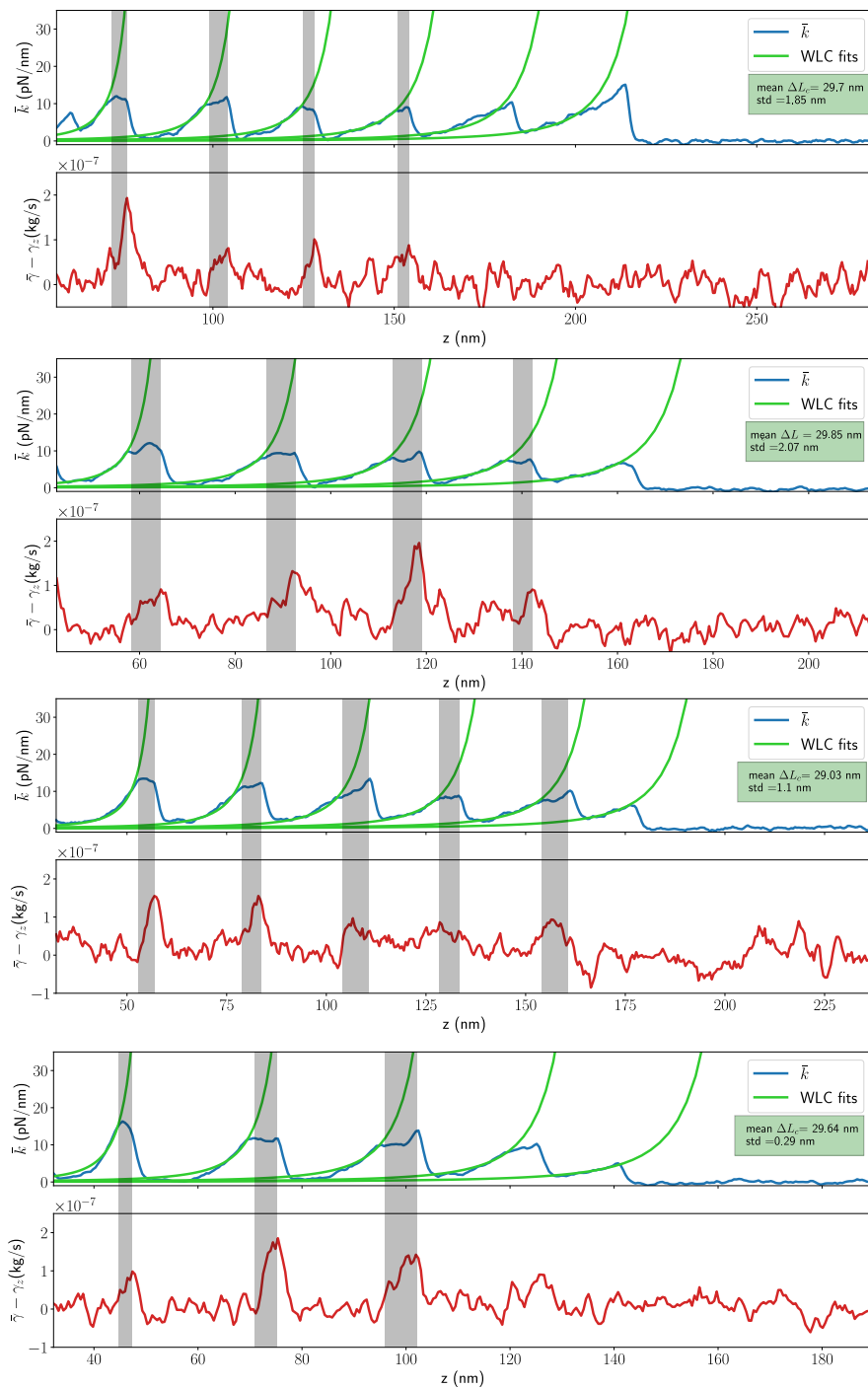


Figure S5:

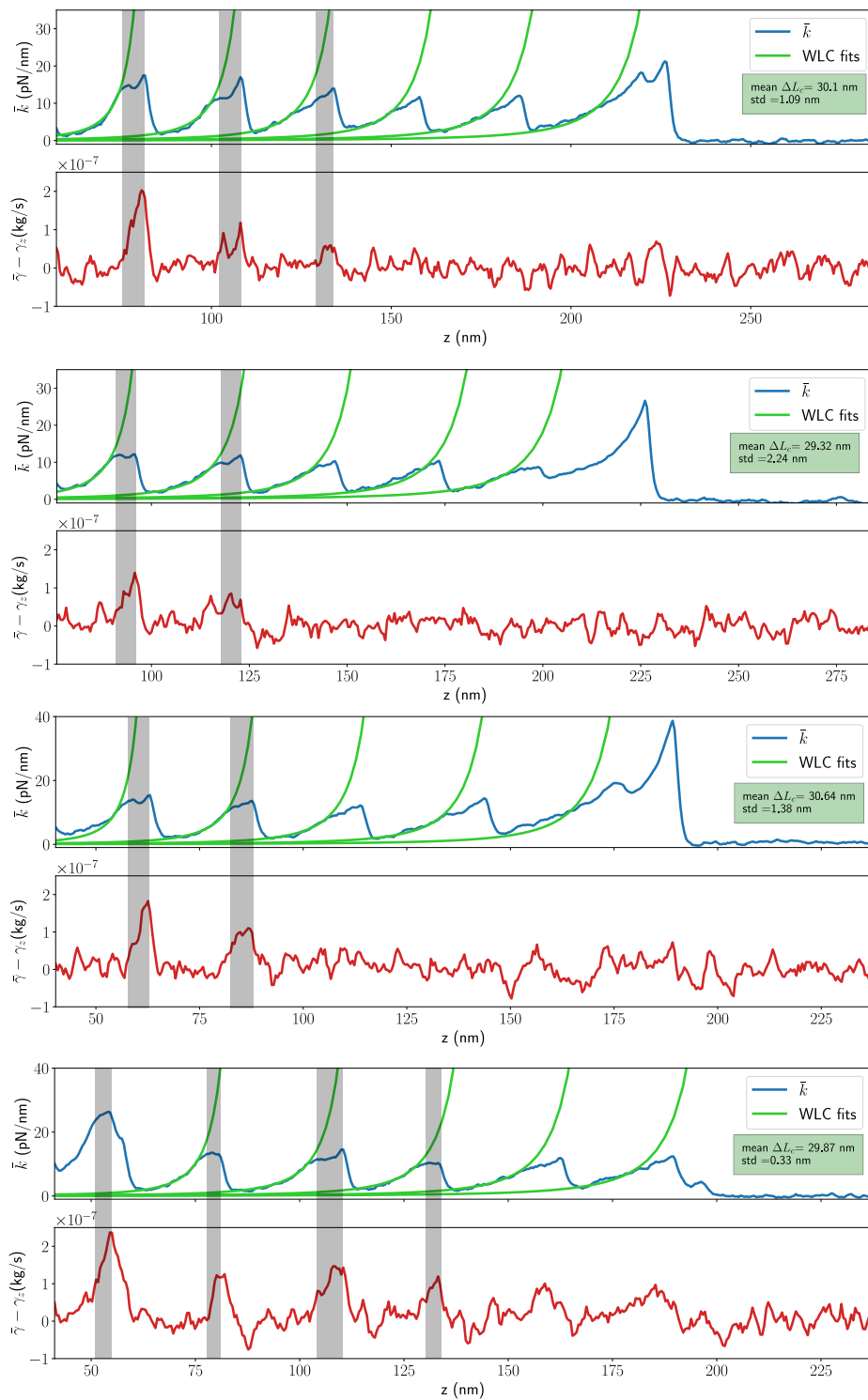


Figure S6:

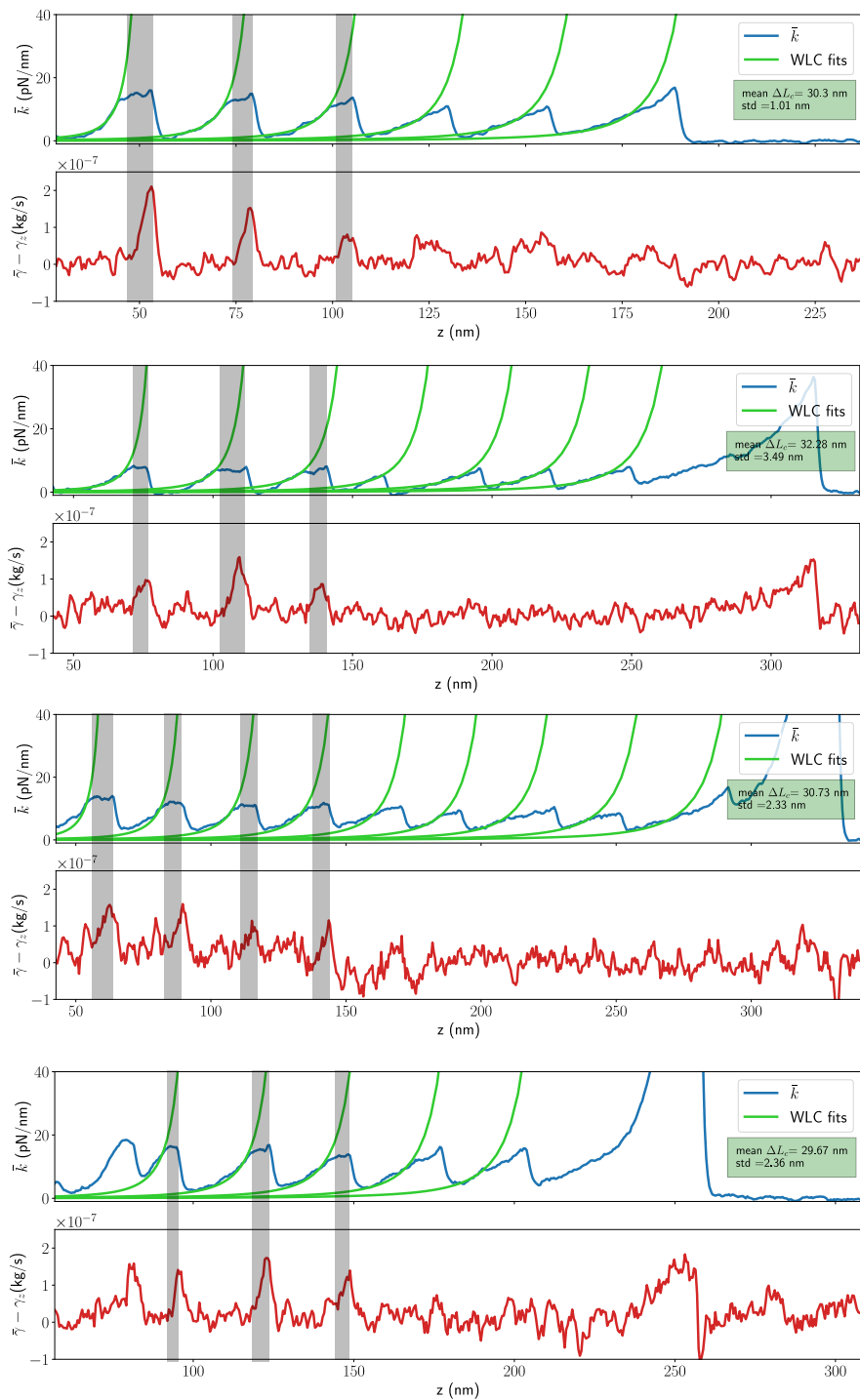


Figure S7: The total stiffness-extension and total friction coefficient-extension profiles. The analysis of these curves yields stiffness and friction coefficient of individual domains, plotted in Fig. 4c and d.

3 Features observed in stiffness-extension and friction coefficient-extension profiles

3.1 The wlc plots of stiffness versus extension and Force versus extension

The force-extension behavior of a polymer chain is described by Marko-Siggia approximation of Worm-Like-Chain (wlc) model.¹⁶ According to the model, which has been shown to work for experimental data obtained using AFM and optical tweezers, the force versus extension is given by

$$F = \frac{k_B T}{p} \left(\frac{1}{4(1 - z/L_c)^2} - \frac{z}{L_c} - \frac{1}{4} \right) \quad (21)$$

The derivative F of this with extension z is stiffness k_{wlc}

$$k_{wlc} = \frac{dF}{dz} = \frac{k_B T}{pL_c} \left(\frac{1}{2(1 - z/L_c)^3} + 1 \right) \quad (22)$$

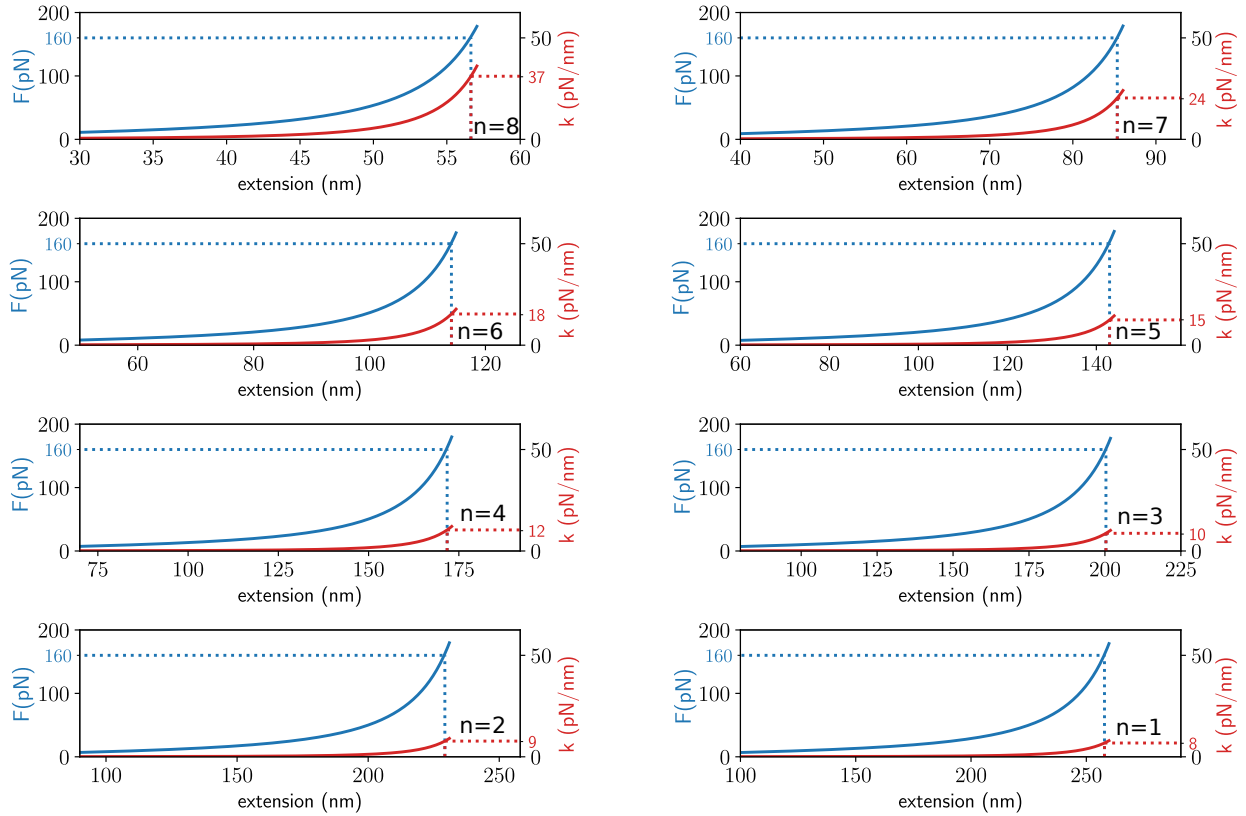


Figure S8: Force and stiffness variation for each contour length. With each unfolding event, the unfolded chain becomes more compliant for the same force (~ 160) pN, thereby dominating the total measured stiffness. The protein is unfolded before the chain has a comparable stiffness to the folded domains so that they start to contribute to the total stiffness measurement. Thus, last four domains unfold without contributing to the stiffness measurement and wlc fits to the entire stretching of the chain.

Fig. S8 shows plots of these two equations for a fixed persistence length $p = 0.4$ nm and different contour lengths L_c . The L_c depends on number of folded domains n . As domains are sequentially unfolded in AFM experiments, each unfolded domain adds 29 nm to the total contour length. The plots clearly indicate that the stiffness of the chain rises sharply compared to the force for larger contour lengths. This means that the chain is more compliant under a given force for larger contour lengths. When there are 8 folded domains, the stiffness is 37 pN/nm for 160 pN, the force at

which the domain unfolds. It is comparable to the domain stiffness. This progressively decreases to 12 pN/nm when there are only 4 domains which are folded. It means that next domain unfolds before the chain becomes stiff enough, so that folded domains contribute to the total stiffness measurement. For last four domains, the contribution to the total stiffness by folded domains is immeasurably low and the wlc model fits to the entire stretching of the molecule till next domain unfolds. It is noteworthy that the dissipation is also close to zero in this region.

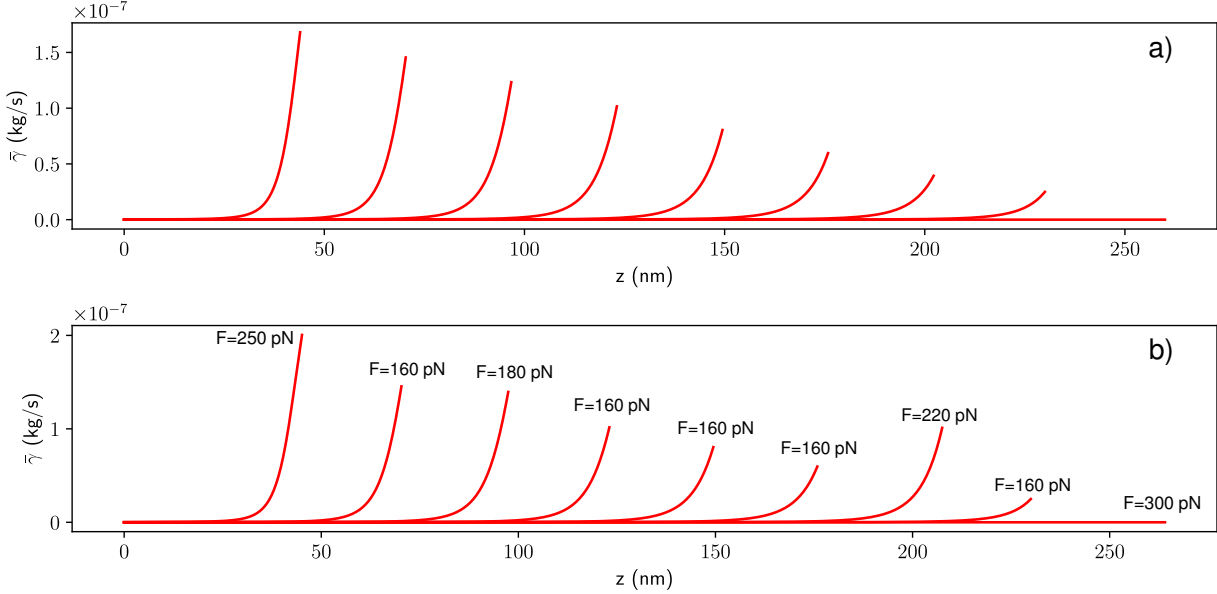


Figure S9: a) The plot of total friction coefficient $\bar{\gamma}$ calculated at different tip-substrate separation z using equation 10. Here the γ_z is taken to be zero. The contour length is increased by 29 nm after each event and n is reduced by one. k_f and γ_f are taken as 100 pN/nm and 2×10^{-6} kg/s. The variation in friction coefficient with respect to z is plotted for each event until 160 pN is reached. This is the average force at which domain unfolds. b) Each domain unfolds at slightly different force as the unfolding events are stochastic. We take similar values for all other parameters as in a, however $\bar{\gamma}$ is now calculated and plotted for each event until unfolding force values are reached for that event. In both plots, it is seen that the peak in the total friction coefficient decreases as domains unfold one by one. The height of the peak also depends on the unfolding force for that event. In experiments, friction coefficient after fourth domain is not observed due to reasons provided in the previous subsection.

3.2 Decreasing peak in total friction coefficient for consecutive unfolding events

It is observed that the total friction coefficient peak height decreases as more and more domains unfold. This feature also needs explanation. It should be noted that the measured damping coefficient is not only of the domain that is unfolding. It is a combined response of all the domains and depends on many other quantities. The $\bar{\gamma}$ given by equation 10, depends not only on friction coefficient of individual domain γ_f , but also on stiffness of the wlc chain k_{wlc} , the number of folded domains n , and stiffness of the individual domain k_f . Of these, we argue that the γ_f and k_f remains fixed. Taking them as constants, we plot equation 10 in Fig. S9a to obtain the total damping coefficient $\bar{\gamma}$ as a function tip-substrate separation till 160 pN is reached and reducing n by 1 after each unfolding event and increasing the L_c by 29 nm. The plot clearly shows that the peak $\bar{\gamma}$ is decreasing with increasing unfolding events. Moreover, we plot the $\bar{\gamma}$ till the unfolding force for the individual unfolding events, if these events are stochastic and the unfolding force for individual domains differ from one another. This results in non-monotonous decrease in the total $\bar{\gamma}$ for each peak as seen in Fig. S9b. The plots explain the observed decrease in the total friction coefficient as the domains sequentially unfold.

4 Fibre-Interferometer based AFM and single molecules

In a typical force spectroscopy experiment, wherein the cantilever is oscillated either at the base or the tip-end, it is treated as a point mass having mass m and spring constant k . Typically this works in ambient or ultra-high vacuum

conditions, where quality factors are large (100 -10000). It has been increasingly realized that such modelling will fail in viscous media such as water and other liquids where the quality factor is low (1-5). When experiments are performed on-resonance and the resultant amplitude and phase are used to quantitatively obtain the stiffness and friction coefficient of the system beneath the tip, it is difficult to predict the phase behaviour seen in the experiments by modelling the cantilever alone.^{31,32} In particular, for single-molecule experiments, it poses additional challenge of controlling the extensions produced in the molecule in a given oscillation cycle.

The above discussion necessitates the off-resonance operation to measure the viscoelasticity of single molecules. There are attempts in the literature to measure single molecule viscoelasticity by driving the cantilever off-resonance.^{30,33-35} However, the cantilever stiffness (and resonance frequency) in all these measurements is low (~ 0.020 pN/nm) and does not allow true off-resonance conditions. In such situations, for a driven cantilever- from either base or tip -, the amplitude reaches its maximum value at resonance and the phase is approximately 90° (Fig. S10). When the tip interacts with a purely elastic element, the resonance of the system shifts to a higher frequency fig.S10b. If one is not operating at *true* off-resonance, such shifts produce a phase signal. If the molecule has internal friction, it is very difficult to separate it from the effects produced due to stiffness change. In experiments especially in a liquid environment, it becomes very difficult to predict how the phase will behave when the interaction contains both conservative and dissipative components.

This issue is resolved if the experiment is performed using stiff cantilevers (0.6 - 1 N/m) having high resonance frequencies (25 - 50 KHz in water). This ensures true off-resonance operation where the phase lag is close to zero and any shift in resonance due to stiffness change does not cause any phase shift. With these cantilevers we get a window of off-resonance condition, in which the measured phase lags are result of damping provided by the molecule alone. See Fig. S10(c) and (d) In single-molecule measurements, the measured stiffness of the poly-protein system is of the order of 10 pN/mm. In true off-resonance, this will produce 1 percent change in tip amplitude. Moreover, it will produce an even smaller bending signal. (For instance with a base amplitude of 1 nm, the tip amplitude is 0.99 nm) This produces a bending amplitude ($A_0 - A$) of 10 pm. Such signals can not be measured using a conventional deflection detection scheme, which is equipped to measure the bending alone . See figure S11 and S12. The interferometer based AFM measures tip displacements. This allows measurement of the extension produced in the molecule directly, in contrast to the deflection detection scheme, which measures bending in the cantilever (Fig S11).

Secondly, a free-of -artefacts, off-resonance operation, is possible with tip excitation and measuring the bending in the cantilever. Here the complications associated with phase response are absent. However, the cantilever and the protein system are in parallel with each other now and the stiffness of the protein system is roughly 100 times smaller than the cantilever (Cantilever: 1 N/m and protein: 10 mN/m). The change in amplitude after protein is attached to the tip is thus 1 part in 100 and is difficult to detect. One needs to use a cantilever having stiffness comparable to that of protein. This makes off-resonance operation difficult as shown in Figure S10. Hence, for measuring the viscoelastic response of protein by exciting the cantilever, interferometer based measurements with base excitation are most suitable. Currently, the measurement bandwidth is limited due to noise floor at low frequencies. See figure S13. This can be overcome by improving laser coherence and filters.³⁶

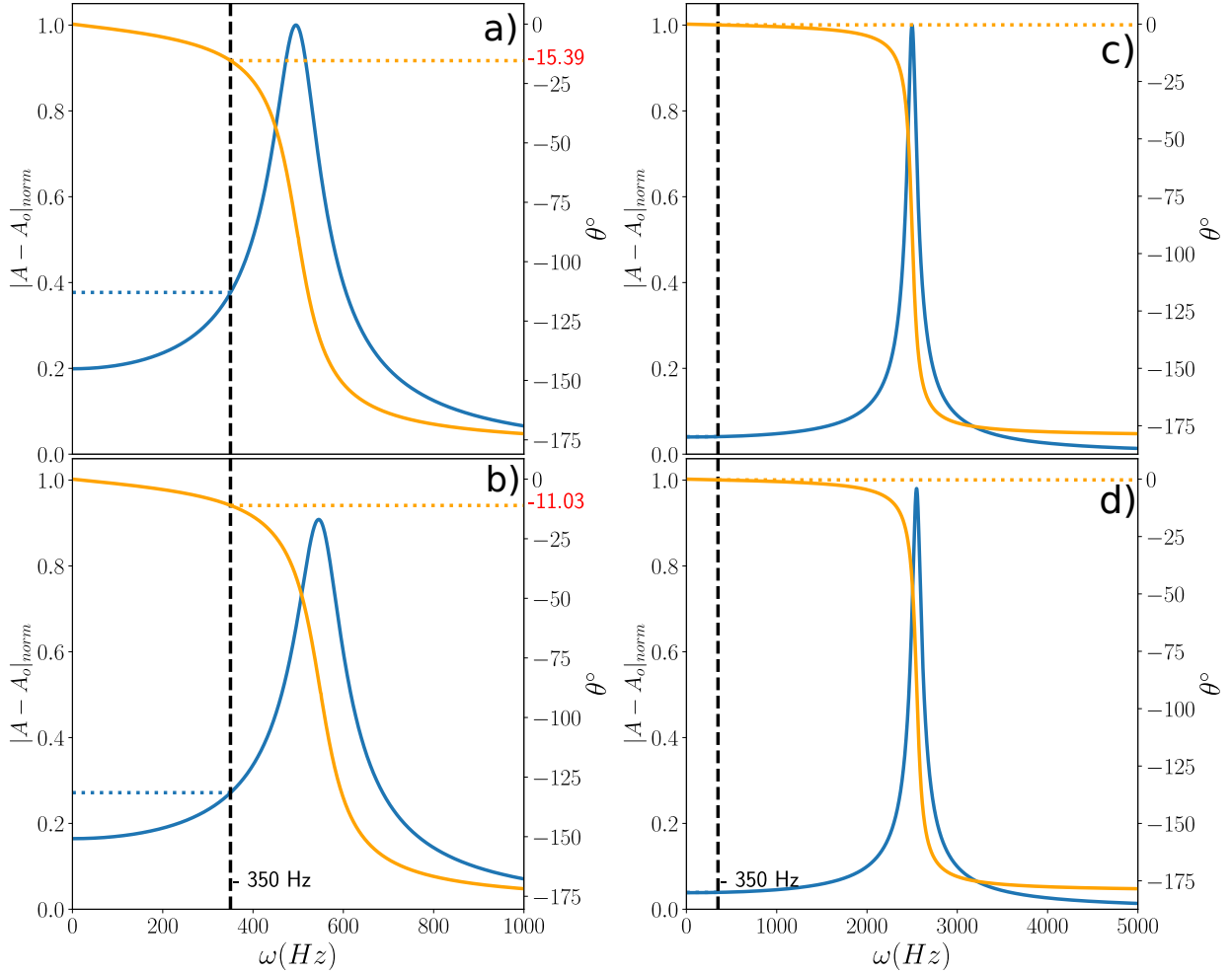


Figure S10: The frequency response of a cantilever which is given by the damped driven harmonic oscillator model: $|A - A_o| = \frac{F_o}{m} \frac{1}{\sqrt{(\omega^2 - \omega_o^2)^2 + (\gamma\omega)^2}}$ and $\tan \theta = -\frac{\gamma\omega}{\omega_o^2 - \omega^2}$ where F_o is the amplitude of the forcing, $|A - A_o|$ represents the bending since A and A_o are the tip and base amplitudes respectively. In the figures $|A - A_o|_{norm}$ is the normalised amplitude. a) and c) have the same parameters except for the resonance frequency, ω_o in c) is 10 times that of a), however the drive frequency is the same for both. In (c) the drive frequency is much lower in comparison to its resonance frequency $\omega_o \gg \omega$ and the amplitude response can be approximated by the static response: $A - A_o \approx \frac{F_o}{k_c}$ and $\theta \approx 0$. In b) and d) an interaction is incorporated that increases the resonance frequency in a) and c) by the same amount. As is evident from the figure, in b) we see a significant change in phase lag, whereas in d) very little change is seen. It is essential to note that one cannot perform off-resonance operation with a cantilever whose response is similar to a), because even at drive frequencies as low as what has been shown here, there is no flat regime in the phase response, any stiffness change can still bring about a change in phase. It is only when stiff cantilevers with high resonance frequencies are used one can achieve this off-resonance operation.

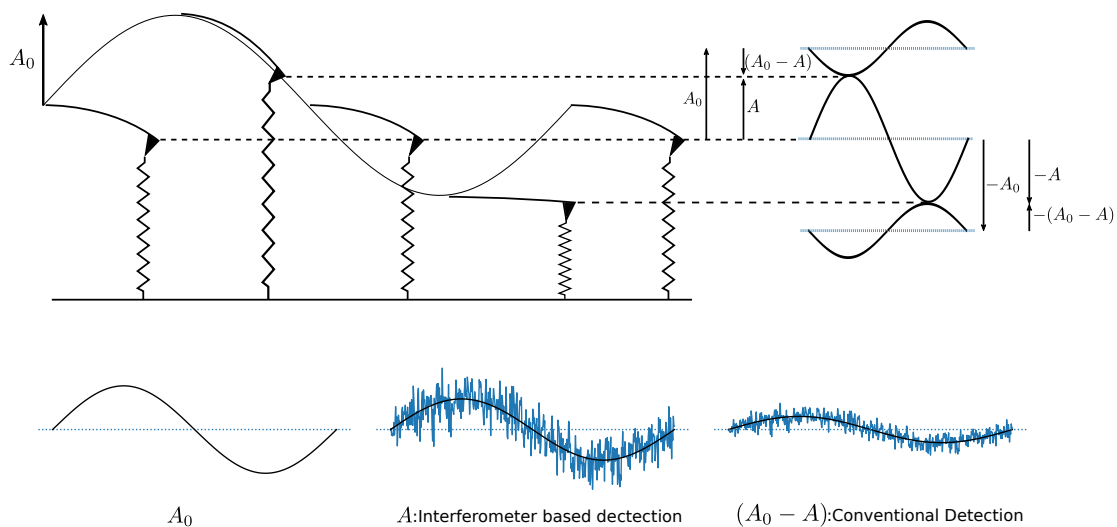


Figure S11: The difference in tip displacement signal (A) and the cantilever bending signal ($A_0 - A$). Interferometer based detection measures tip displacements whereas the deflection detection scheme measures the bending. For stiff cantilevers (0.5 to 1 N/m), the bending is extremely small and difficult to detect. See fig. S12 for experimental data which shows the actual bending and displacement signals. Furthermore, A also represents the amplitude of extension in the molecule over the oscillation, which is directly measured by the interferometer.

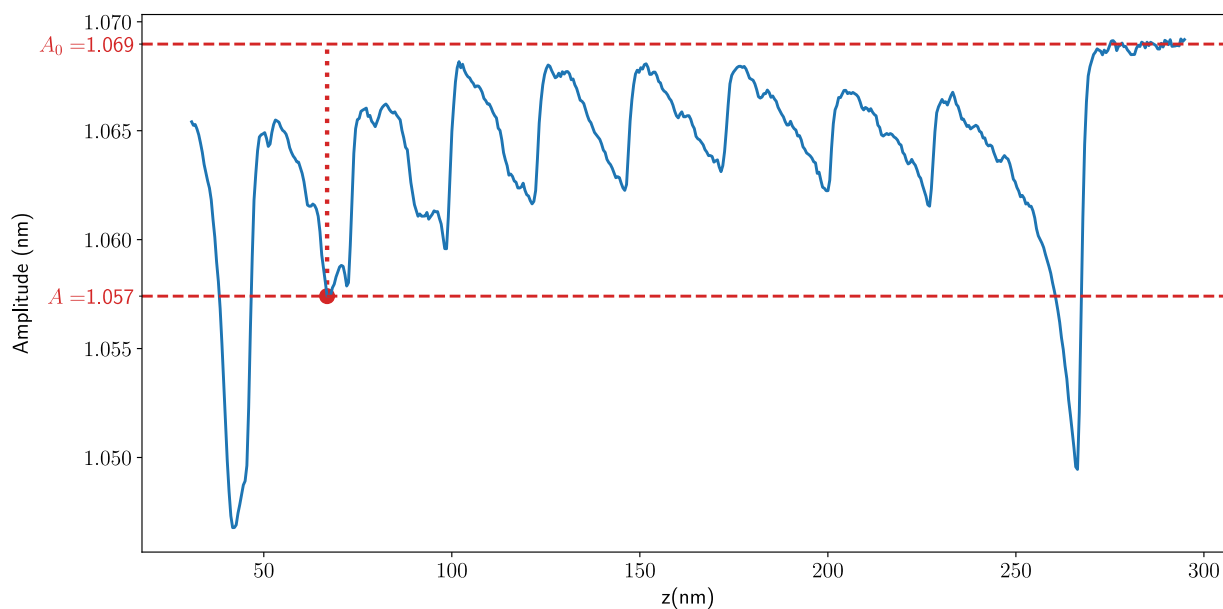


Figure S12: The amplitude A change measured as the cantilever is pulled away from the surface with the protein attached. A_0 is determined from the measured amplitude when the protein detaches from the cantilever. The bending signal in this measurement at the point marked in red $A_0 - A$ is 12 pm. The sensitivity of the interferometer for this measurement was 300 mV/nm.

5 Measurement at lower frequency

In order to see the effect of performing measurements at lower frequencies, the experiments are carried out at 433 Hz. The stiffness data once again shows that at high stretch (shaded region in Fig.S8a) there is contribution to stiffness from folded domains. The stiffness-extension data begins to deviate from the wlc description. We use our analysis in this region to separate the stiffness of folded intermediates from the combined total stiffness. The stiffness of folded intermediates is similar to measurement performed at ~ 2 KHz. The friction coefficient signal, on the other hand is featureless. However the noise floor at this frequency is much higher compared to the frequency window around 2 KHz. It can not be concluded from this data if the internal friction coefficient depends on operational frequency.

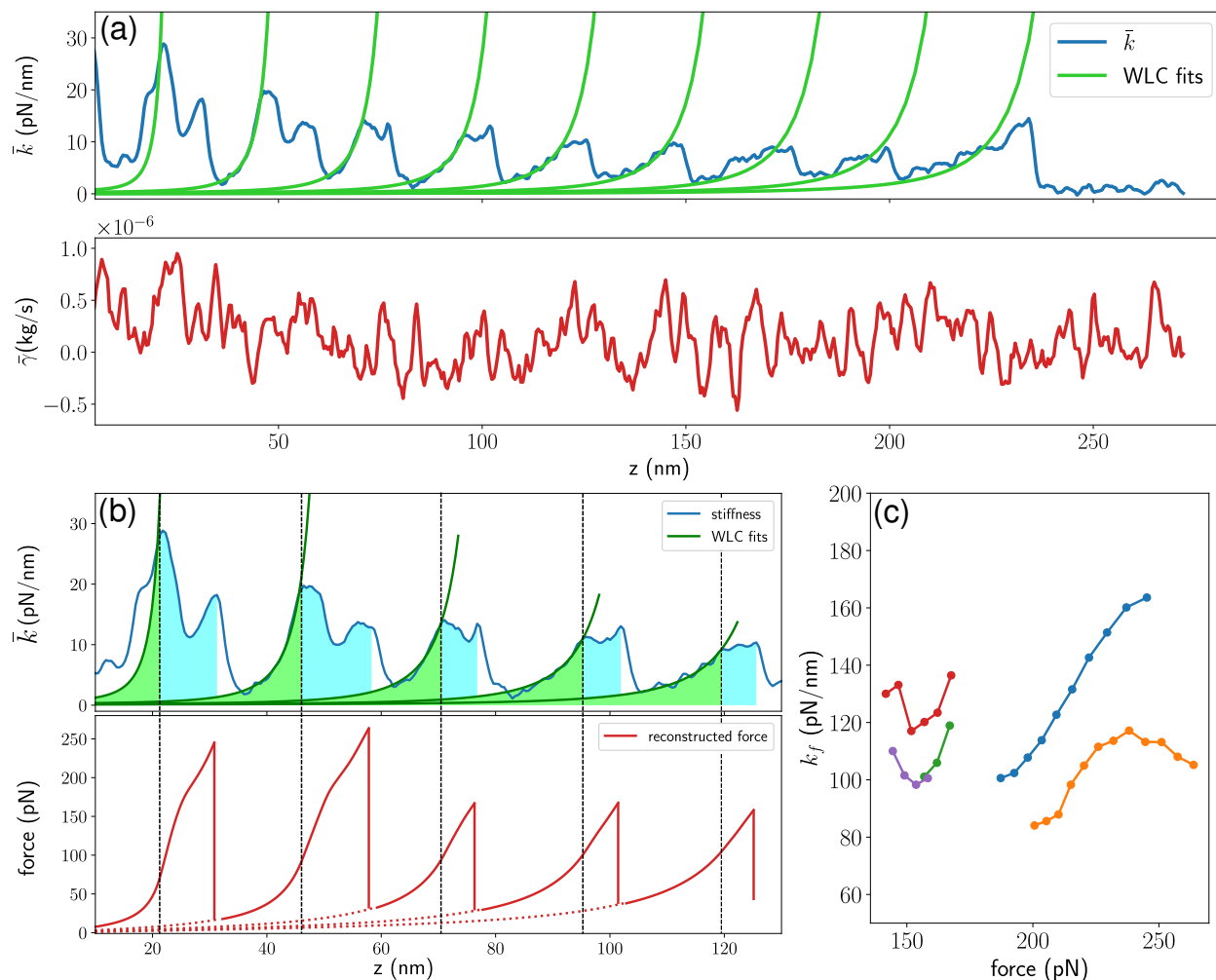


Figure S13: Measured stiffness and friction coefficient for a measurement performed at 433Hz. The noise floor is higher than the measured friction coefficient at 2.1 kHz. The dissipation is hence not measured. The analysis shows that the force profiles and its intersection with dotted vertical lines once again indicates that the transition to region, where stiffness from folded intermediates contribute to the measurement occurs at ~ 95 pN. The stiffness of the intermediate at different forces are also shown. The stiffness values are similar to the experiment performed at higher frequencies around 2 KHz.